



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Bupivacaine
---------------	-------------

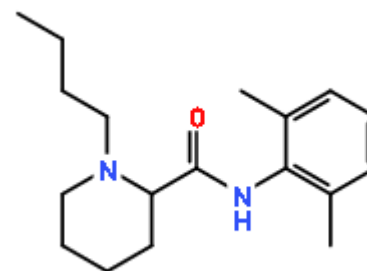

CHROMATOGRAM

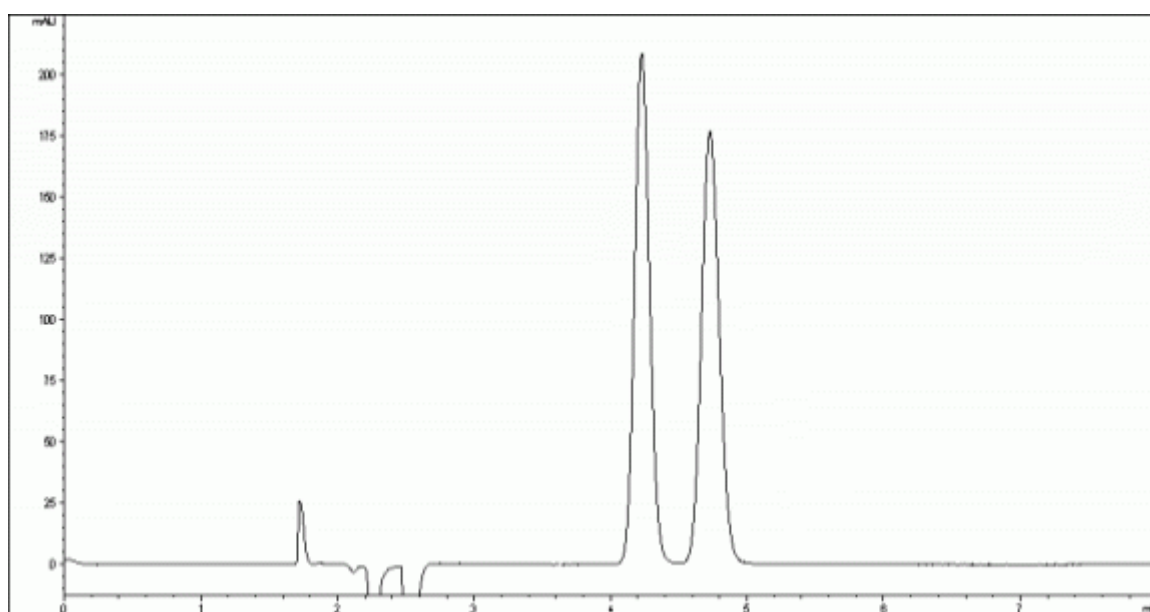

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALCEL® OZ-3 150 x 4.6mm / 3µm
Mobile Phase Composition	n-hexane / 2-propanol / diethylamine = 90 / 10 / 0.1
Chromatographic Mode	Normal Phase Coated CSP

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 230 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	4.23	1.5	1.2	2.26
2	4.73	1.8		