


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Hexobarbital

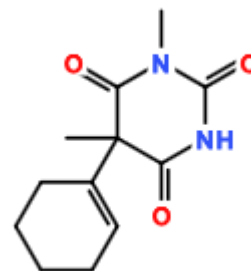

CHROMATOGRAM

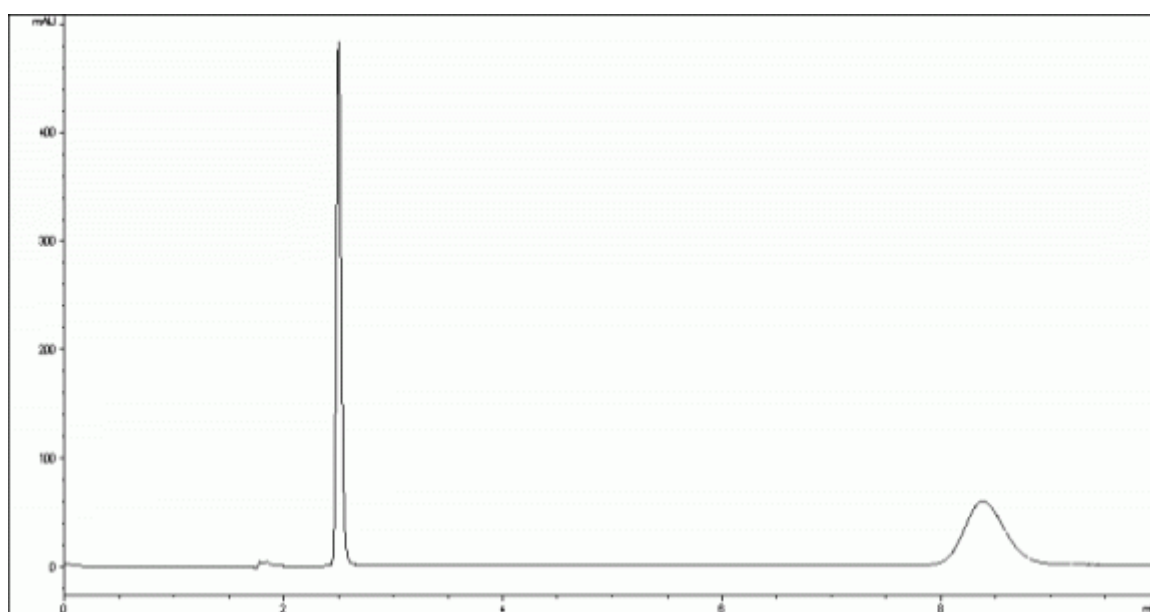

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AY-3 150 x 4.6mm / 3µm
Mobile Phase Composition	n-hexane / ethanol = 50 / 50
Chromatographic Mode	Normal Phase Coated CSP

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 230 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	2.5	0.41	1 - 2	9.02
2	8.38	3.74		15.49