



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	alpha-Methoxyphenylacetic acid
---------------	--------------------------------

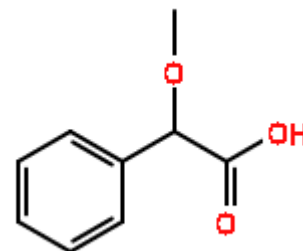

CHROMATOGRAM

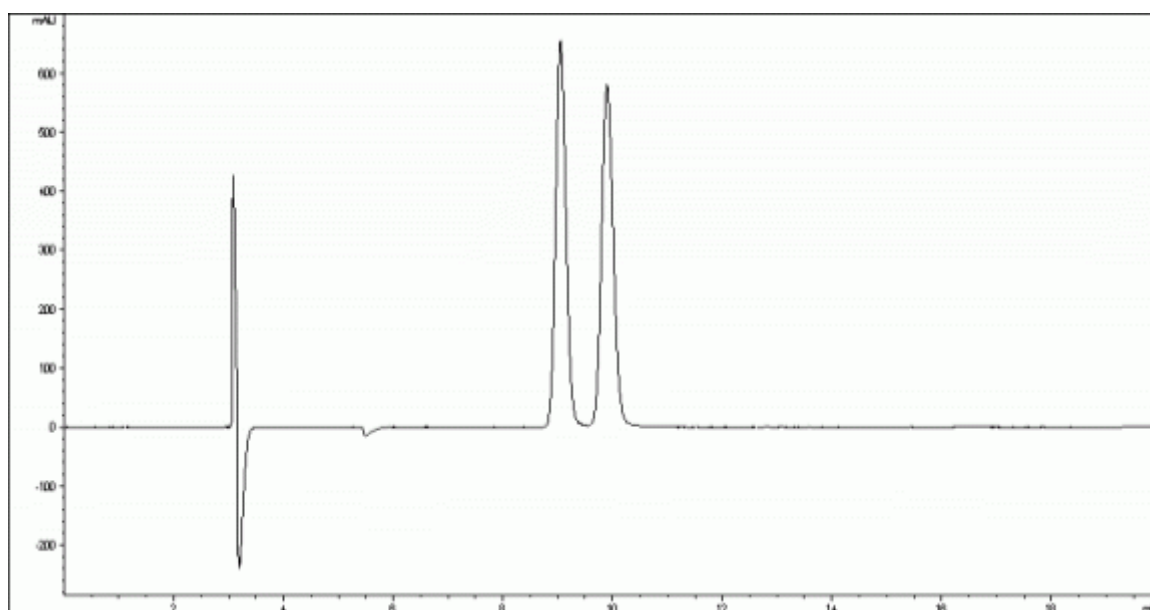

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AZ-H 250 x 4.6mm / 5µm
Mobile Phase Composition	n-hexane / 2-propanol / trifluoroacetic acid = 95 / 5 / 0.1
Chromatographic Mode	Normal Phase Coated CSP

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 210 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	9.05	2.05	1 - 2	2.34
2	9.9	2.34	1.14	