



## APPLICATION RESULTS

## COMPOUND IDENTIFICATION

Compound Name	2-Chloromandelic acid
---------------	-----------------------

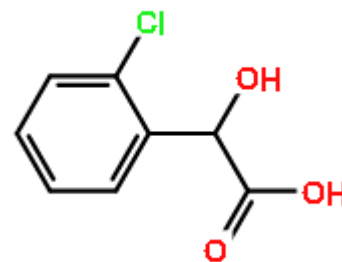

## CHROMATOGRAM

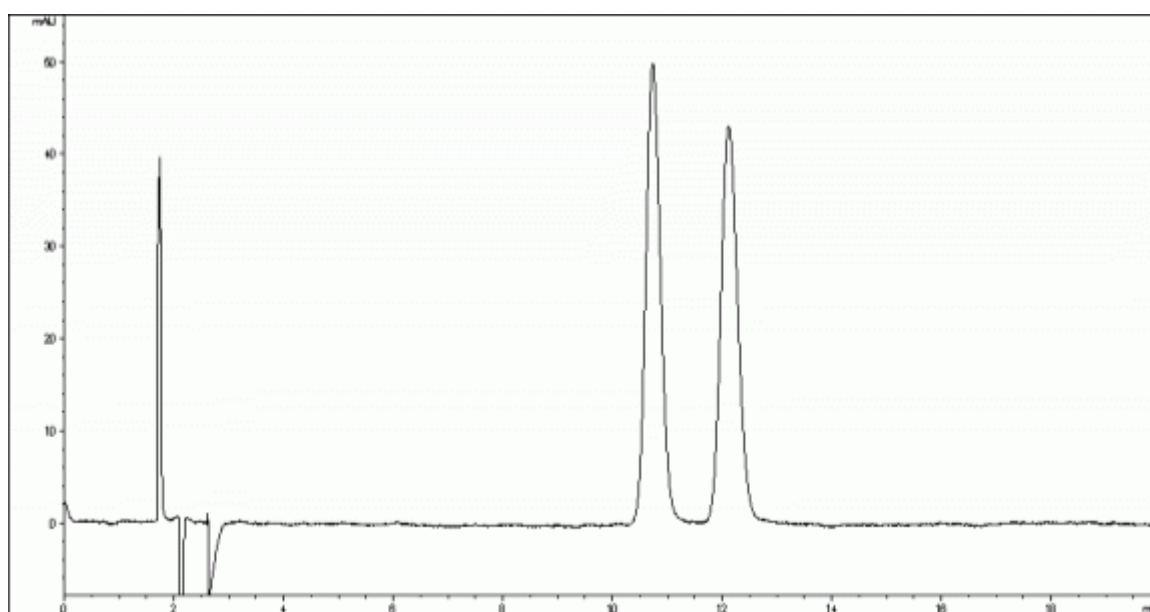

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AS-3 150 x 4.6mm / 3µm
Mobile Phase Composition	n-hexane / 2-propanol / trifluoroacetic acid = 90 / 10 / 0.1
Chromatographic Mode	Normal Phase Coated CSP

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 230 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	10.74	5.39	1 - 2	2.62
2	12.12	6.21		