



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	4-Benzoyloxy-2-azetidinone
---------------	----------------------------

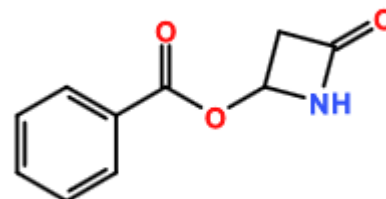

CHROMATOGRAM

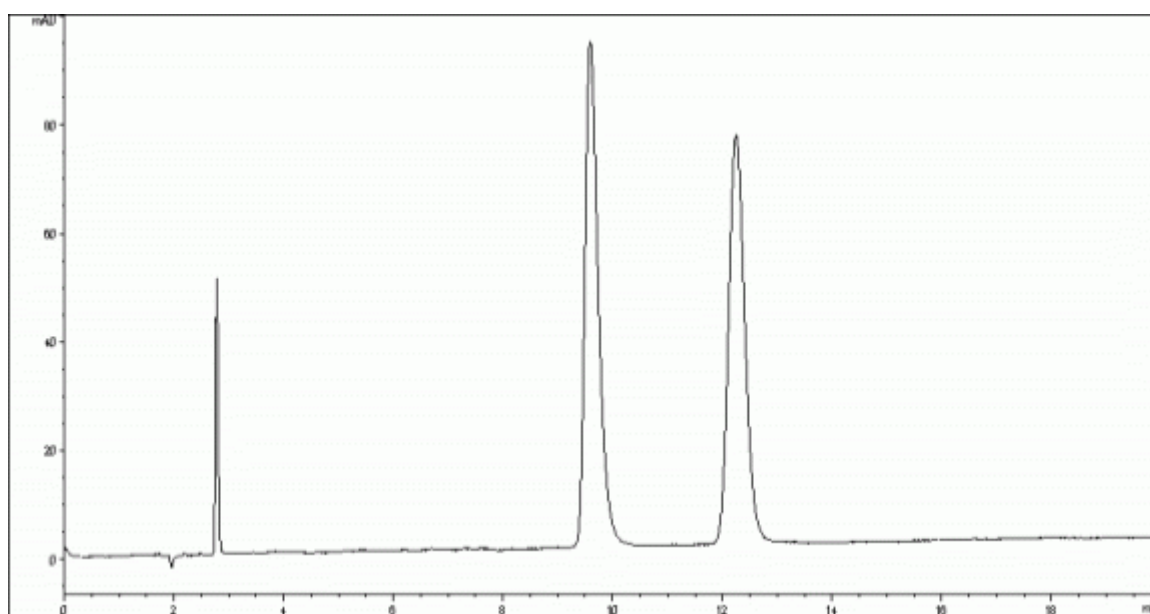

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALCEL® OD-3 150 x 4.6mm / 3µm
Mobile Phase Composition	n-hexane / 2-propanol = 90 / 10
Chromatographic Mode	Normal Phase Coated CSP

Flow rate (ml/min)	1	Injection amount (mg)	0.001
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 230 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	9.59	4.84	1 - 2	5.5
2	12.25	6.46		