


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Metalaxyl

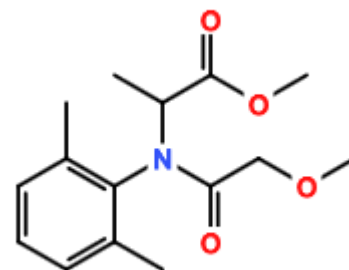

## CHROMATOGRAM

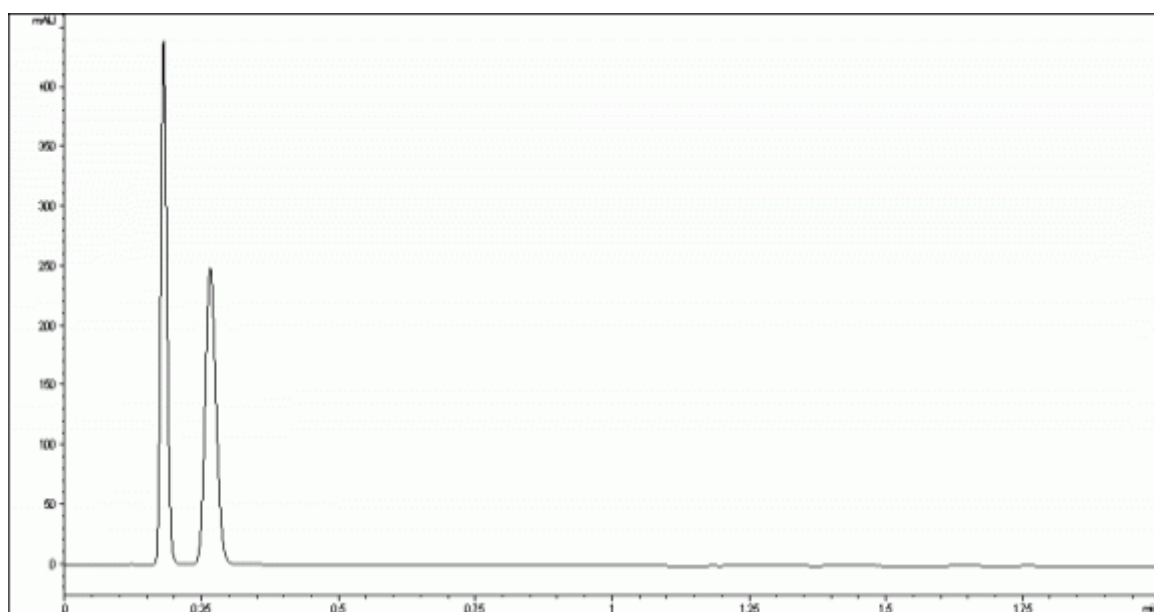

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALCEL® OD-3 50 x 4.6mm / 3µm
Mobile Phase Composition	n-hexane / ethanol = 70 / 30
Chromatographic Mode	Normal Phase Coated CSP

Flow rate (ml/min)	5	Injection amount (mg)	0.001
Temperature (°C)	40	Pressure (MPa)	
Detection	UV 230 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	0.18	0.62	1 - 2	2.23
2	0.27	1.38		