



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Pentobarbital sodium
---------------	----------------------

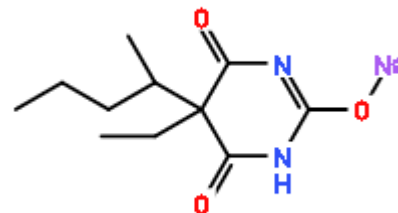

CHROMATOGRAM

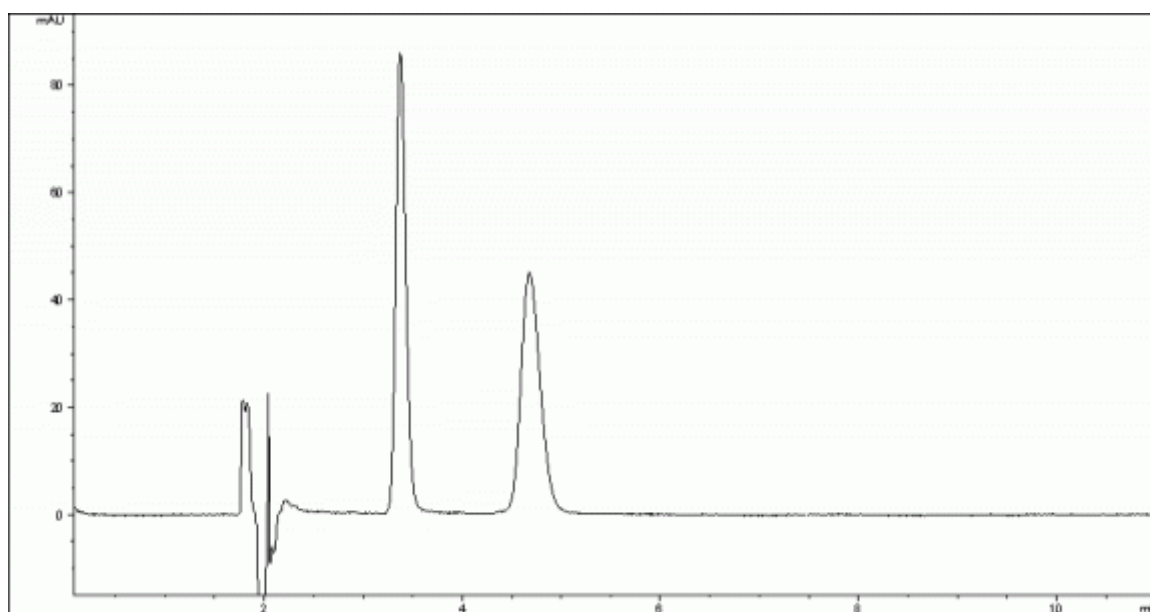

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALCEL® OJ-3 150 x 4.6mm / 3µm
Mobile Phase Composition	n-hexane / 2-propanol / trifluoroacetic acid = 80 / 20 / 0.1
Chromatographic Mode	Normal Phase Coated CSP

Flow rate (ml/min)	1	Injection amount (mg)	0.01
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 230 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	3.37	0.88	1 - 2	4.86
2	4.68	1.61		