



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	trans-Stilbene Oxide
---------------	----------------------

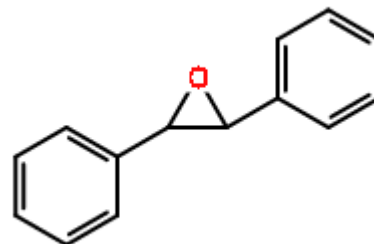

CHROMATOGRAM

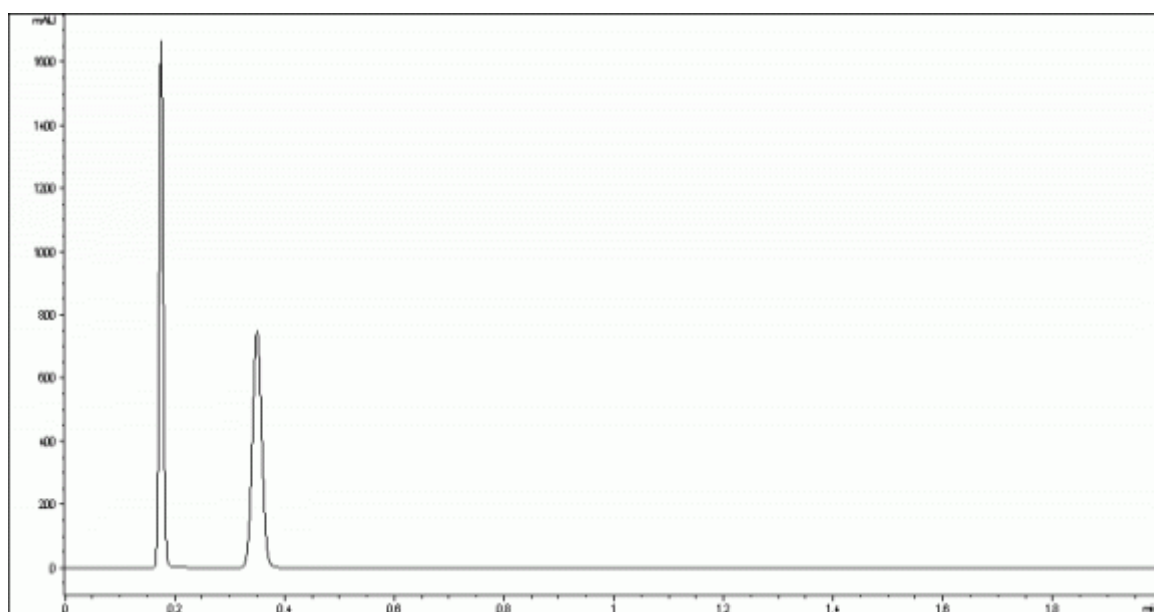

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AD-3 50 x 4.6mm / 3µm
Mobile Phase Composition	n-hexane / 2-propanol = 90 / 10
Chromatographic Mode	Normal Phase Coated CSP

Flow rate (ml/min)	6	Injection amount (mg)	0.002
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 230 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	0.18	0.88	1 - 2	7.72
2	0.35	2.76		