



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	(±)-Hydrobenzoin
---------------	------------------

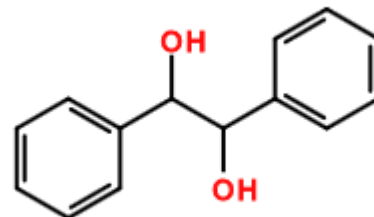

CHROMATOGRAM

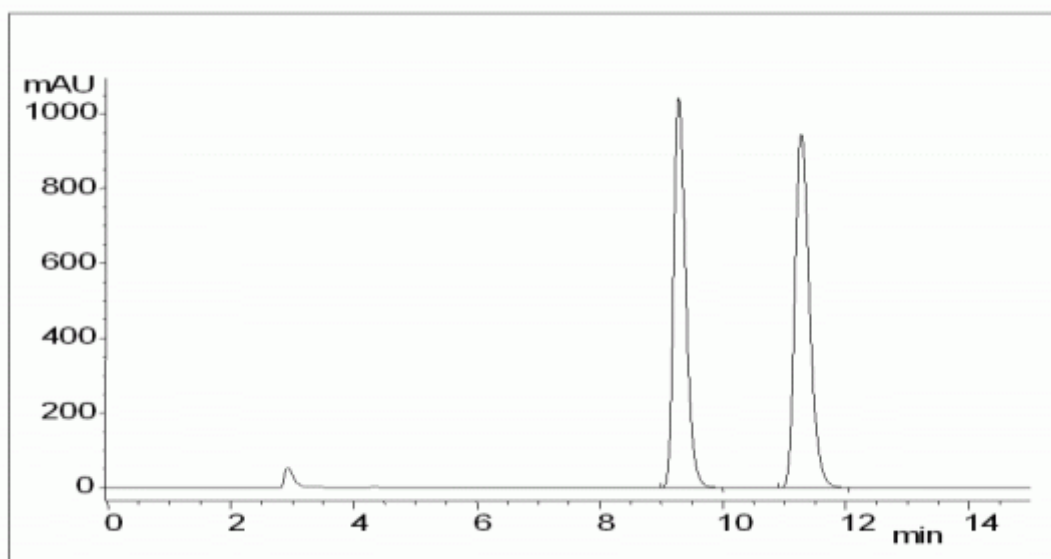

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AD-H 250 x 4.6mm / 5µm
Mobile Phase Composition	n-hexane / ethanol = 85 / 15
Chromatographic Mode	Normal Phase Coated CSP

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 220 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	9.28	2.09	1 - 2	5.11
2	11.28	2.76		