



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	4-Amino-2-chlorodiphenylacetonitrile
---------------	--------------------------------------

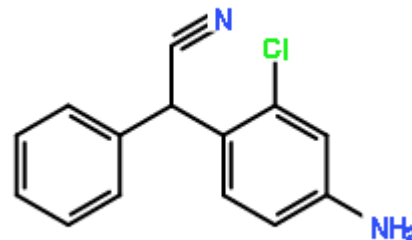

CHROMATOGRAM

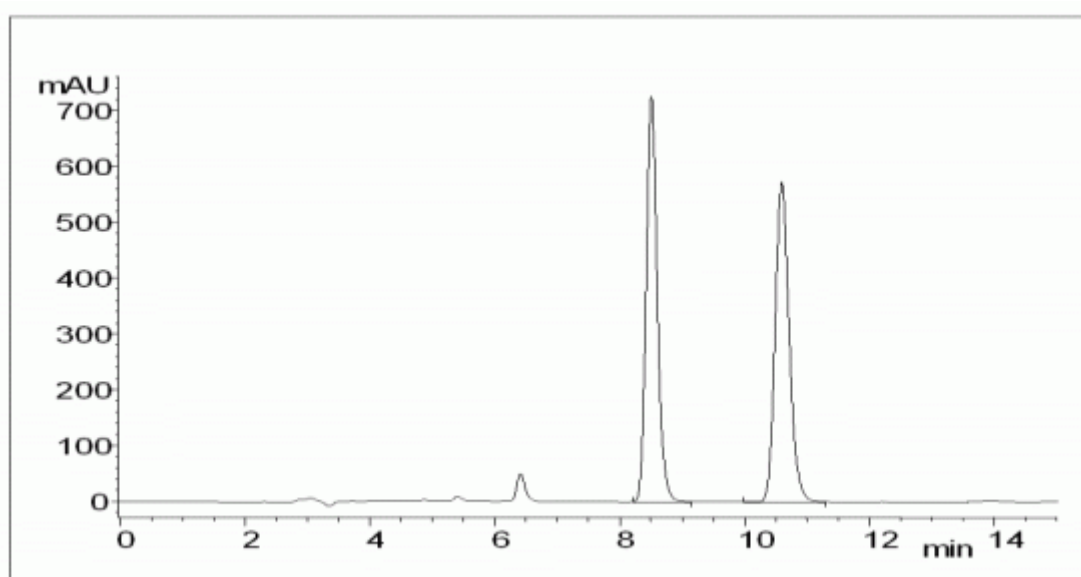

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AD-H 250 x 4.6mm / 5µm
Mobile Phase Composition	n-hexane / 2-propanol = 70 / 30
Chromatographic Mode	Normal Phase Coated CSP

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 210 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	8.51	1.84	1 - 2	5.91
2	10.6	2.53		