



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	4-Amino-2-chlorodiphenylacetonitrile
---------------	--------------------------------------

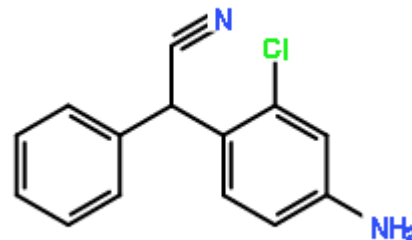

CHROMATOGRAM

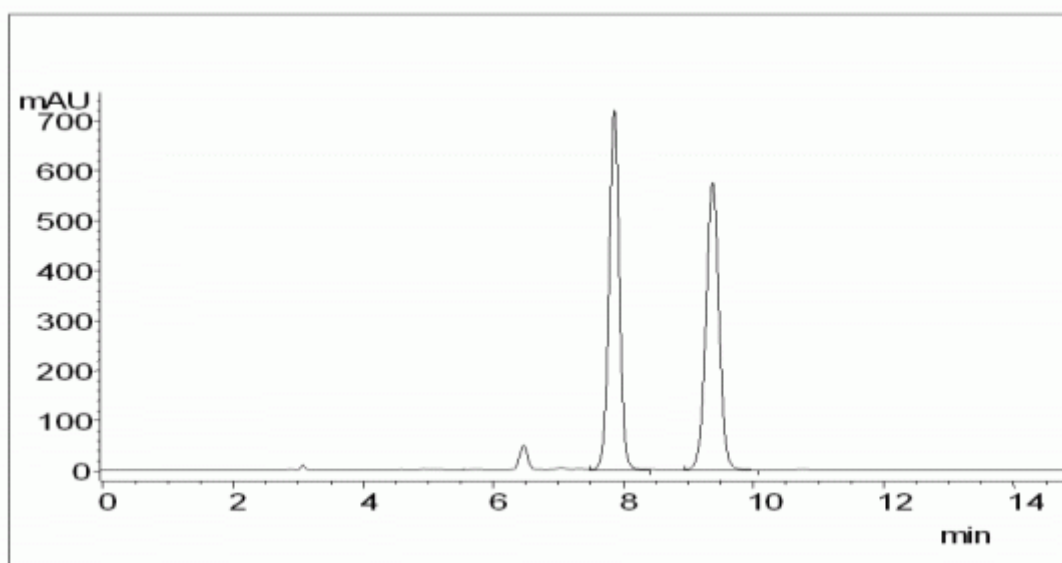

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AS-H 250 x 4.6mm / 5µm		
Mobile Phase Composition	n-hexane / ethanol = 75 / 25		
Chromatographic Mode	Normal Phase Coated CSP		

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 210 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	7.86	1.62	1 - 2	4.68
2	9.38	2.13		