



## APPLICATION RESULTS

## COMPOUND IDENTIFICATION

Compound Name	4-Benzoyloxy-2-azetidinone
---------------	----------------------------

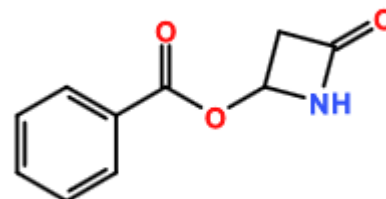

## CHROMATOGRAM

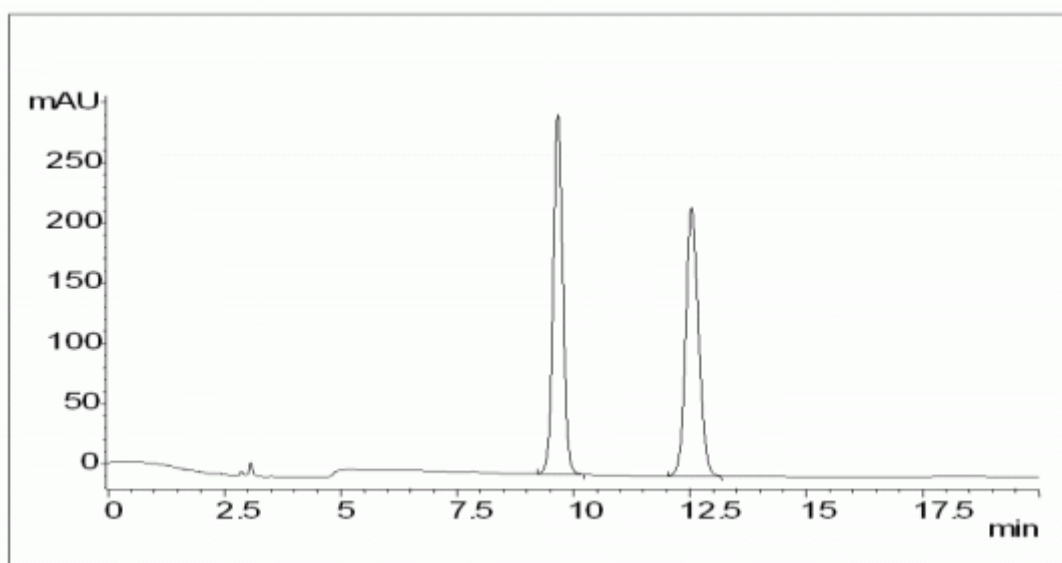

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AS-H 250 x 4.6mm / 5µm
Mobile Phase Composition	n-hexane / ethanol = 75 / 25
Chromatographic Mode	Normal Phase Coated CSP

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 210 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	9.67	2.22	1 - 2	6.6
2	12.54	3.18		