



## APPLICATION RESULTS

## COMPOUND IDENTIFICATION

Compound Name	Alprenolol
---------------	------------

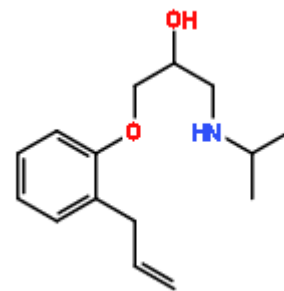

## CHROMATOGRAM

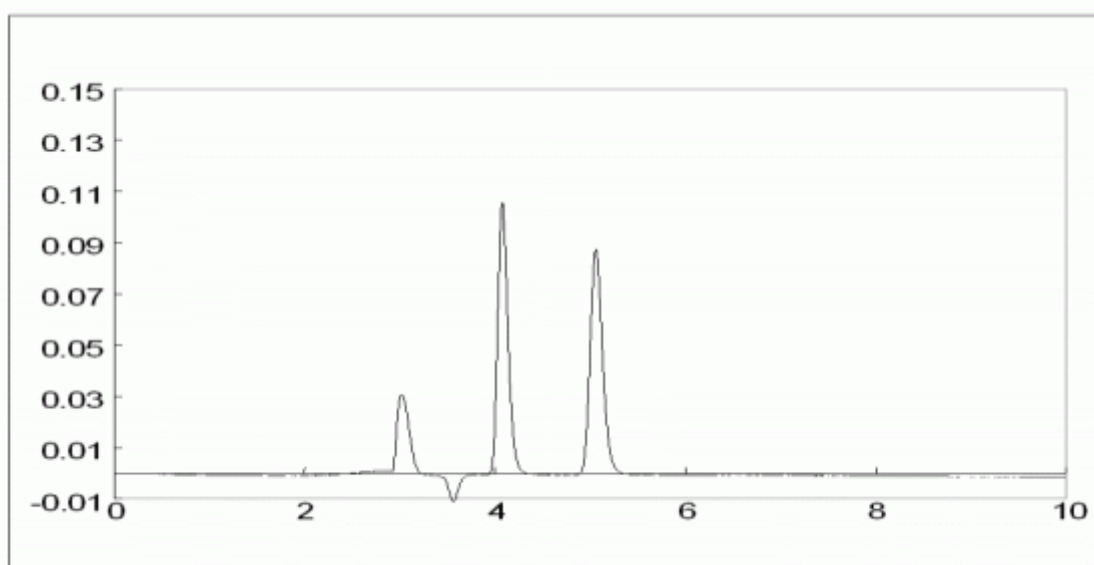

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALCEL® OD-H 250 x 4.6mm / 5µm
Mobile Phase Composition	n-heptane / ethanol / diethylamine = 90 / 10 / 0.1
Chromatographic Mode	Normal Phase Coated CSP

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 254 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	4.07	0.35	1 - 2 1.93	4.52
2	5.05	0.67		