



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Carprofen
---------------	-----------

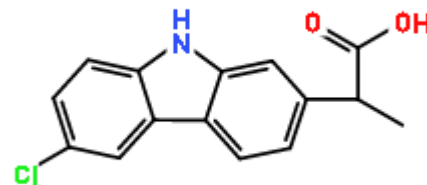

CHROMATOGRAM

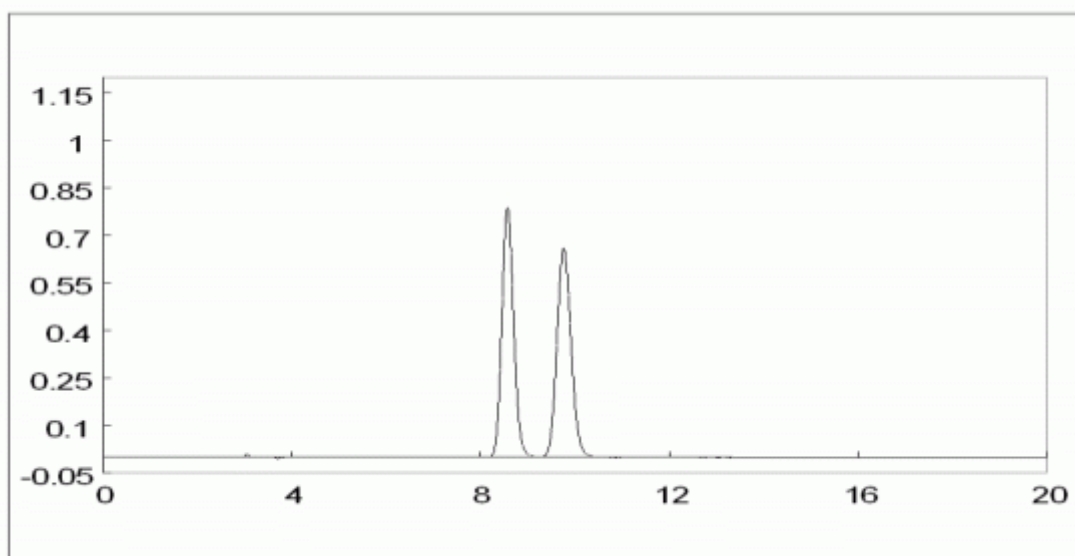

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALCEL® OD-H 250 x 4.6mm / 5µm
Mobile Phase Composition	n-heptane / ethanol / trifluoroacetic acid = 90 / 10 / 0.1
Chromatographic Mode	Normal Phase Coated CSP

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 254 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	8.58	1.86	1 - 2	1.21
2	9.77	2.26		2.35