


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Citalopram

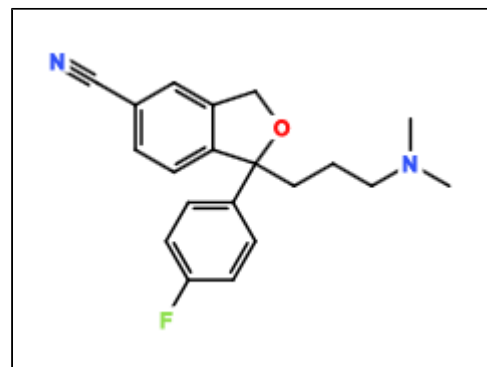

CHROMATOGRAM

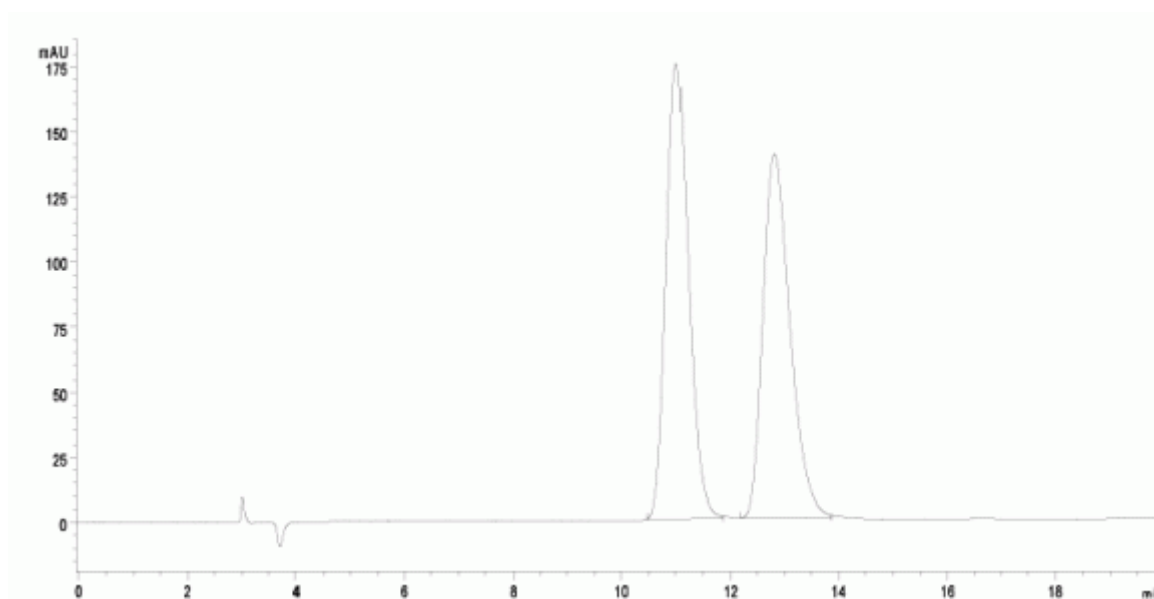

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALCEL® OD-H 250 x 4.6mm / 5µm
Mobile Phase Composition	n-heptane / 2-propanol / diethylamine = 94 / 6 / 0.1
Chromatographic Mode	Normal Phase Coated CSP

Flow rate (ml/min)	1	Injection amount (mg)	0.01
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 250 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	11	2.66	1 - 2	1.23
2	12.8	3.27		2.2