


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Hexobarbital

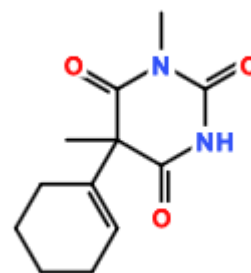

CHROMATOGRAM

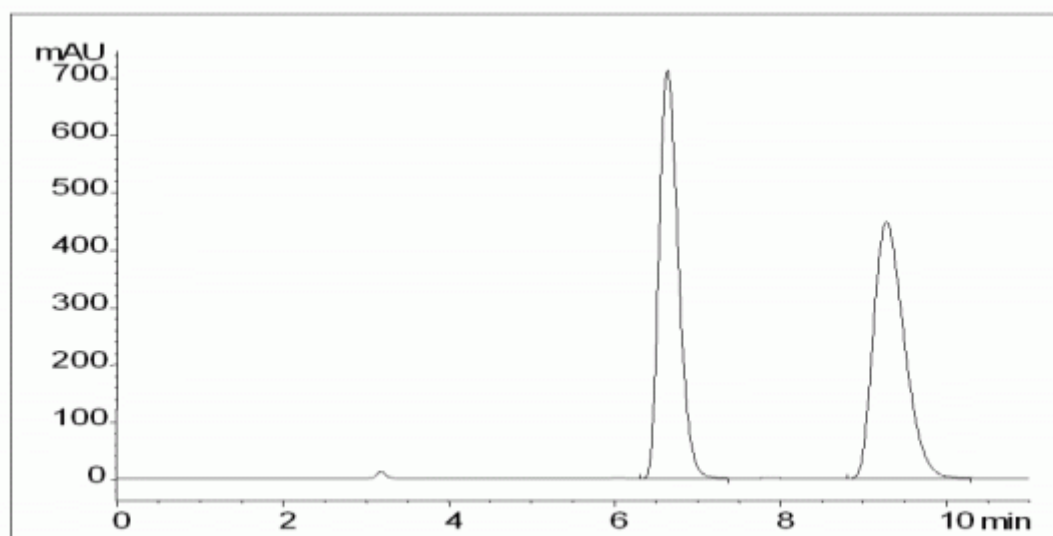

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AS-H 250 x 4.6mm / 5µm
Mobile Phase Composition	n-heptane / 2-propanol = 80 / 20
Chromatographic Mode	Normal Phase Coated CSP

Flow rate (ml/min)	1	Injection amount (mg)	0.01
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 210 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	6.63	1.21	1 - 2	4.7
2	9.27	2.09		