


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

DL-Leu-DL-Leu

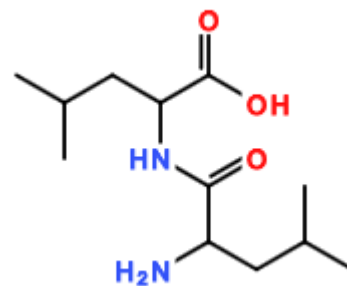

## CHROMATOGRAM

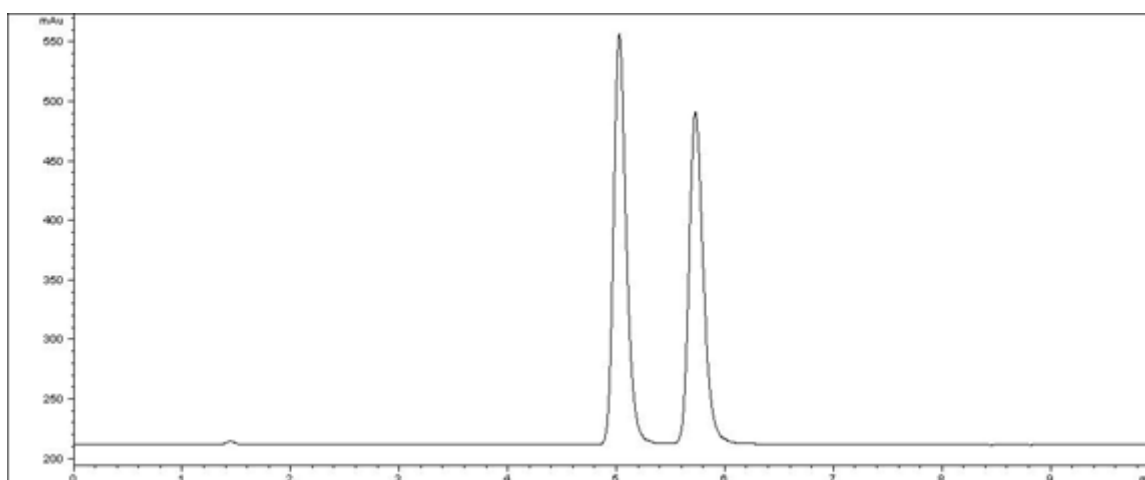

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® ZWIX(+) 150 x 3mm / 3µm
Mobile Phase Composition	50mM formic acid + 25mM diethylamine in methanol / tetrahydrofuran / water = 49 / 49 / 2
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	0.0025
Temperature (°C)	25	Pressure (MPa)	8.7
Detection	ELSD	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	5.02	1.81	1 - 2	3.19
2	5.72	2.2		