


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Mephesisin

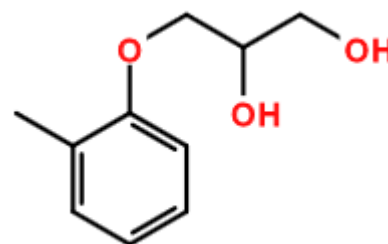

CHROMATOGRAM

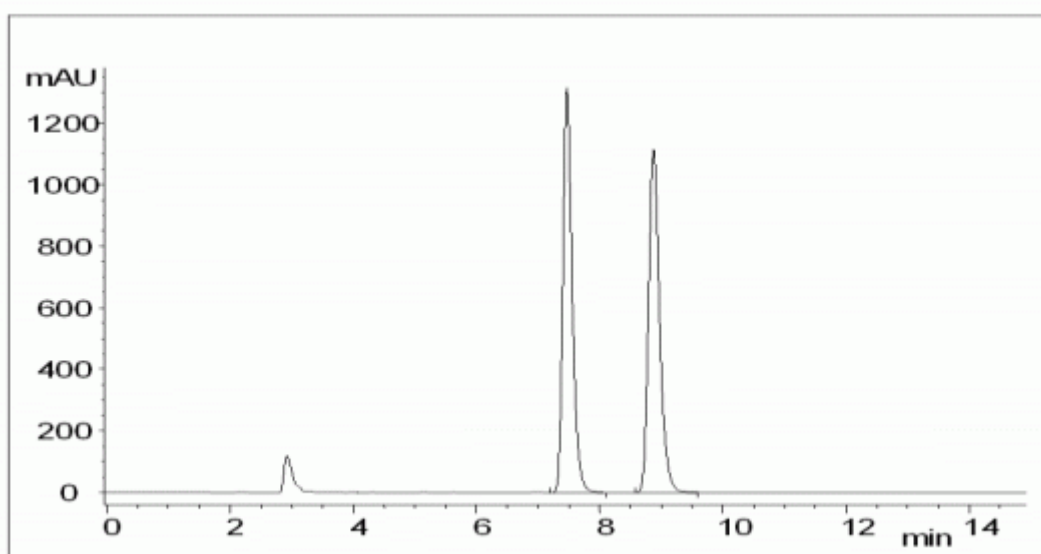

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AD-H 250 x 4.6mm / 5µm
Mobile Phase Composition	n-hexane / ethanol = 85 / 15
Chromatographic Mode	Normal Phase Coated CSP

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 210 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	7.47	1.49	1 - 2	4.82
2	8.88	1.96		