



## APPLICATION RESULTS

## COMPOUND IDENTIFICATION

Compound Name	Nicardipine
---------------	-------------

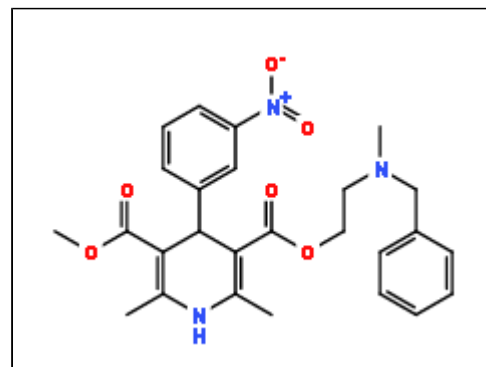

## CHROMATOGRAM

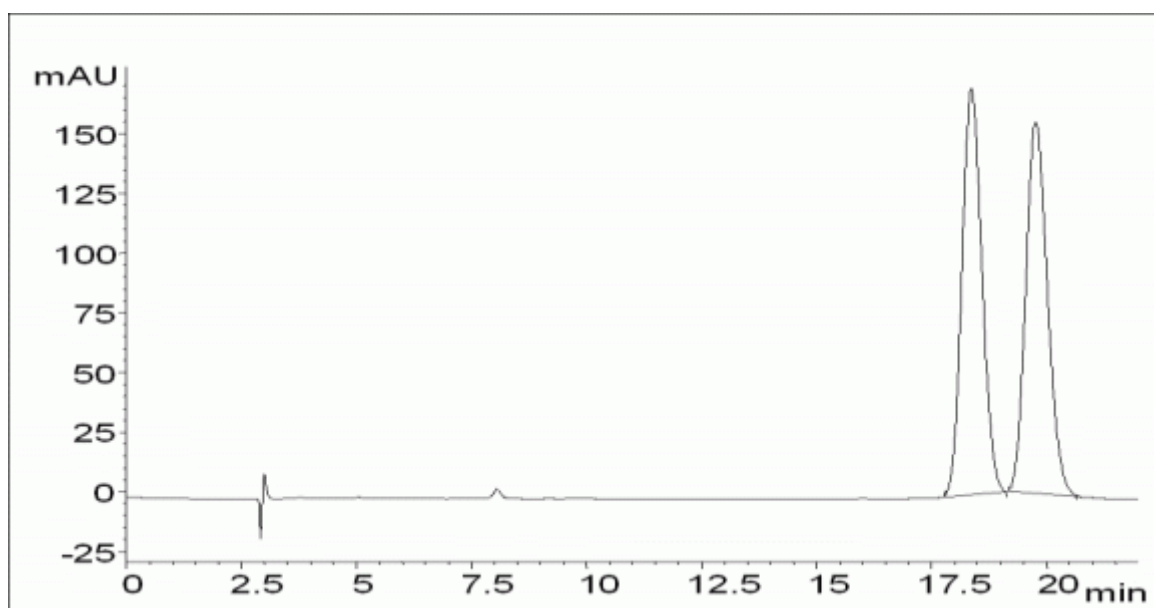

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AD-H 250 x 4.6mm / 5µm
Mobile Phase Composition	n-hexane / 2-propanol / diethylamine = 93 / 7 / 0.1
Chromatographic Mode	Normal Phase Coated CSP

Flow rate (ml/min)	1	Injection amount (mg)	0.01
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 275 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	18.36	5.12	1 - 2	1.1
2	19.76	5.58		1.7