


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Glutethimide

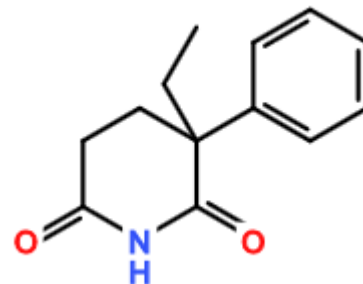

## CHROMATOGRAM

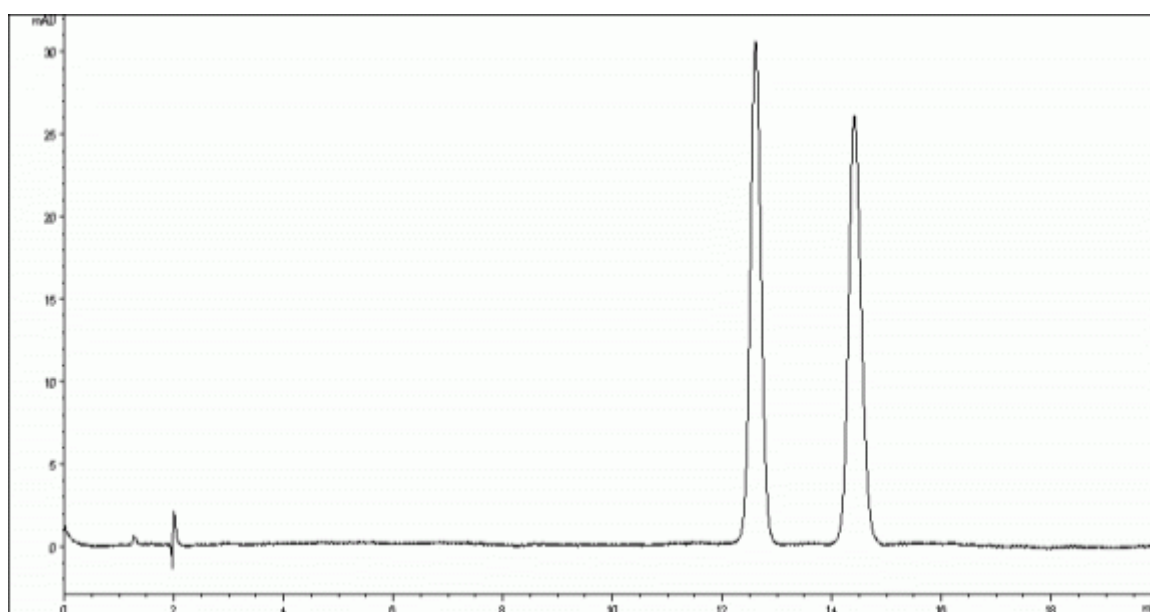

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IB-3 150 x 4.6mm / 3µm
Mobile Phase Composition	water / acetonitrile = 70 / 30
Chromatographic Mode	Reverse Phase

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 230 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	12.61	6.2	1 - 2	4.62
2	14.41	7.24		