



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	DL-Lysine
---------------	-----------

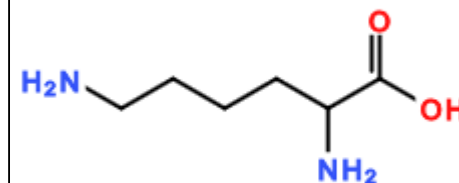

CHROMATOGRAM

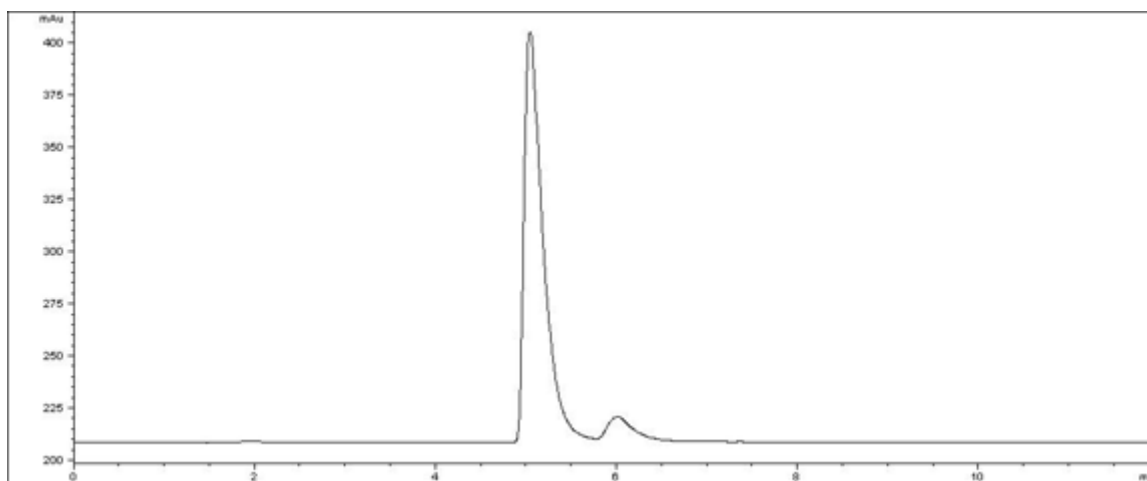

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® ZWIX(-) 150 x 3mm / 3µm
Mobile Phase Composition	50mM formic acid + 25mM diethylamine in methanol / water = 98 / 2
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	0.0025
Temperature (°C)	25	Pressure (MPa)	
Detection	ELSD	Elution order	D/L

No.	Rt	K'	Alpha Value	Resolution
1	5.04	1.8	1 - 2	2.19
2	6	2.33		