



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	N-3,5-DNB-DL-Leucine
---------------	----------------------

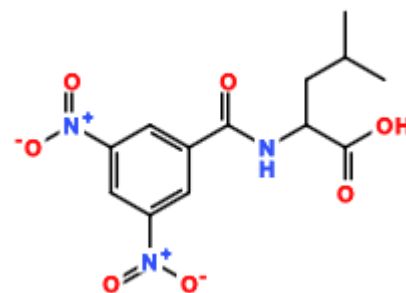

CHROMATOGRAM

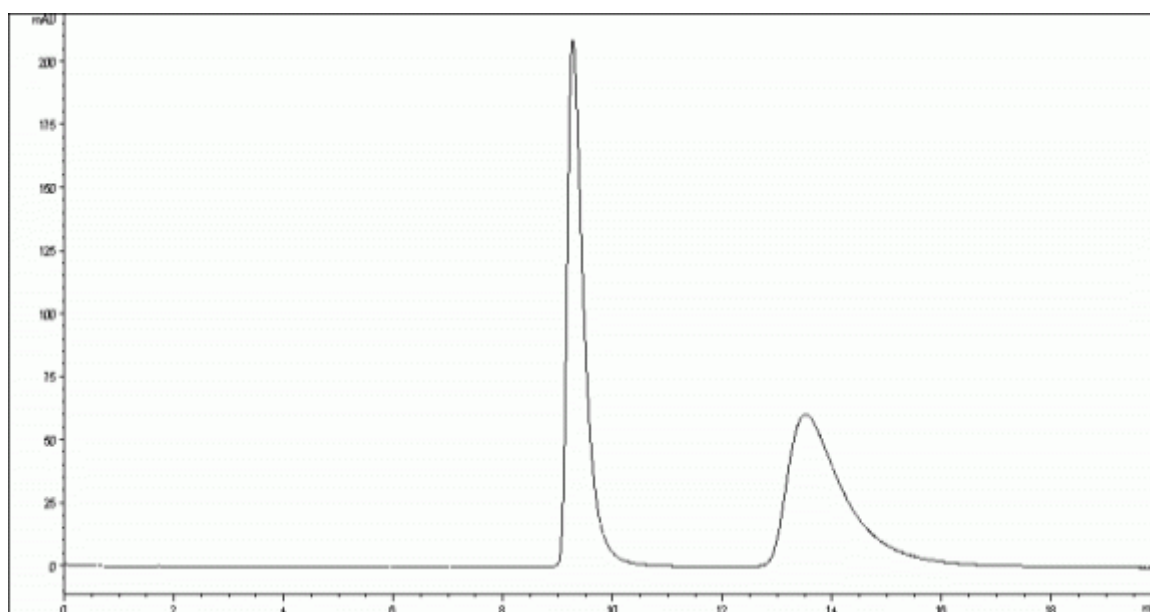

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IA-3 150 x 4.6mm / 3µm
Mobile Phase Composition	trifluoroacetic acid aq. pH2.0 / acetonitrile = 60 / 40
Chromatographic Mode	Reverse Phase

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 254 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	9.27	4.3	1 - 2	3.75
2	13.53	6.74		