


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Glutethimide

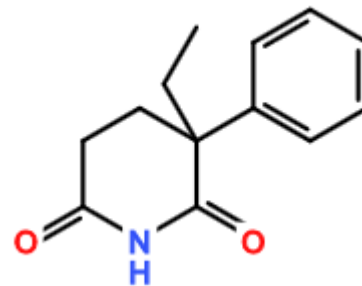

## CHROMATOGRAM

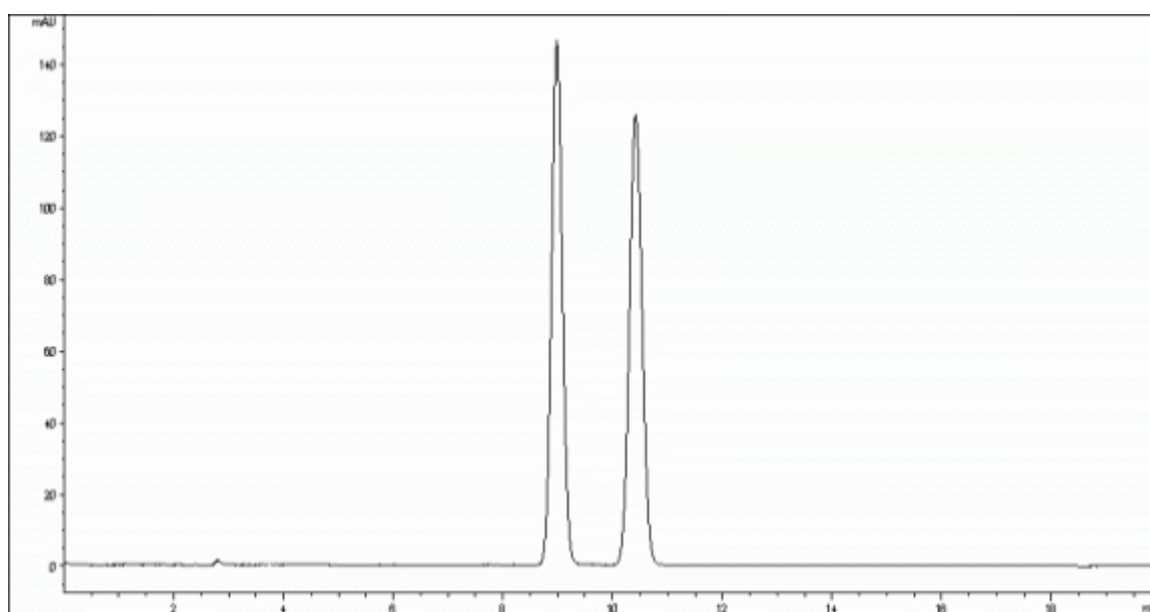

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALCEL® OD-3R 150 x 4.6mm / 3µm
Mobile Phase Composition	water / acetonitrile = 60 / 40
Chromatographic Mode	Reverse Phase

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 220 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	8.99	4.52	1.2	3.65
2	10.42	5.4		