



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Naproxen
---------------	----------

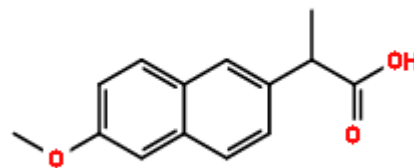

CHROMATOGRAM

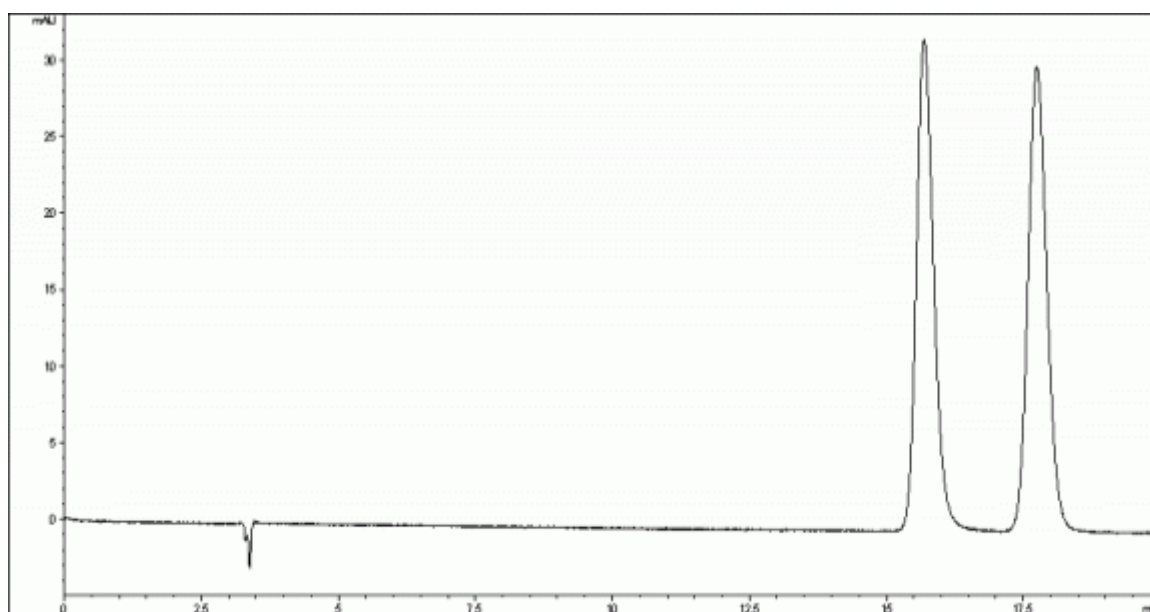

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AD-3R 150 x 4.6mm / 3µm
Mobile Phase Composition	formic acid aq. pH2.0 / acetonitrile = 60 / 40
Chromatographic Mode	Reverse Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	0.001
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 254 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	15.69	3.51	1 - 2	3.47
2	17.75	4.1		