


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Benzoin

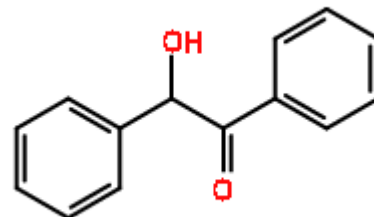

CHROMATOGRAM

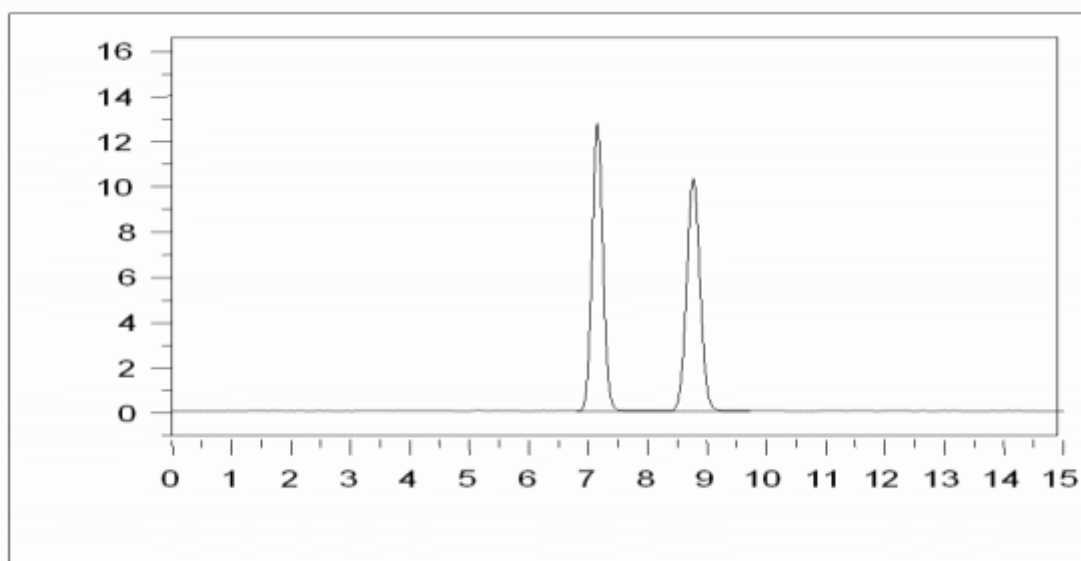

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AS-RH 150 x 4.6mm / 5µm
Mobile Phase Composition	water / acetonitrile = 60 / 40
Chromatographic Mode	Reverse Phase

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 254 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	7.16	3.09	1 - 2	1.3
2	8.77	4.01		