



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Pinocebrin-7-methyl ether
---------------	---------------------------

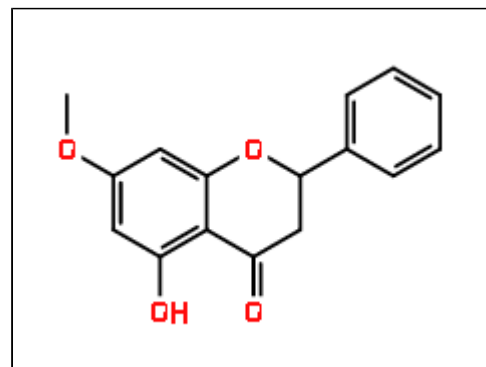

CHROMATOGRAM

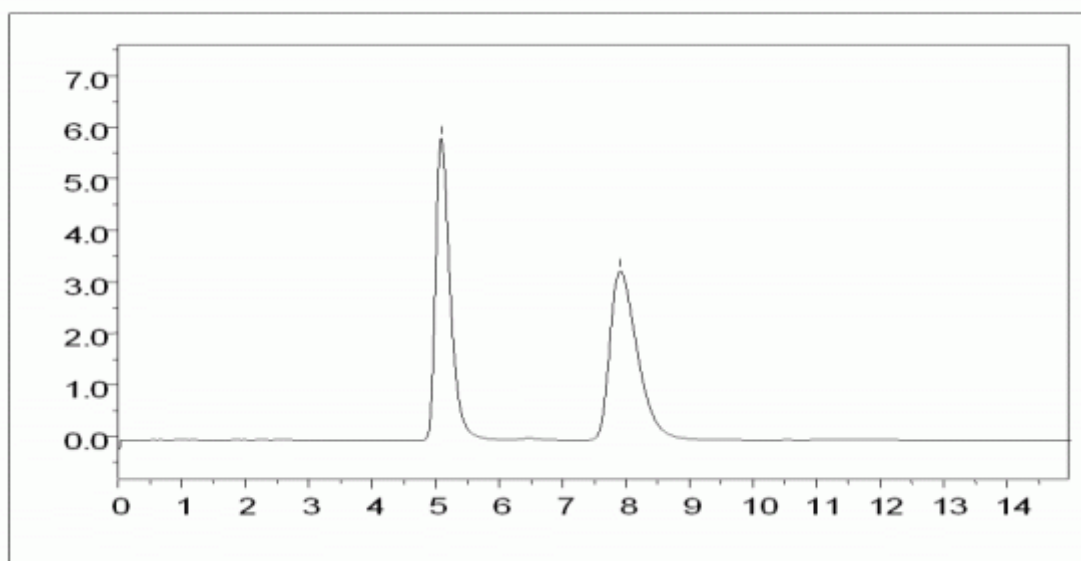

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AD-RH 150 x 4.6mm / 5µm
Mobile Phase Composition	water / acetonitrile = 20 / 80
Chromatographic Mode	Reverse Phase

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 254 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	5.09	3.04	1 - 2	4.74
2	7.91	5.28		