



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Temazepam
---------------	-----------

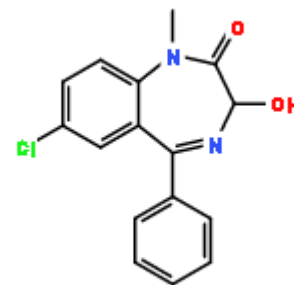

CHROMATOGRAM

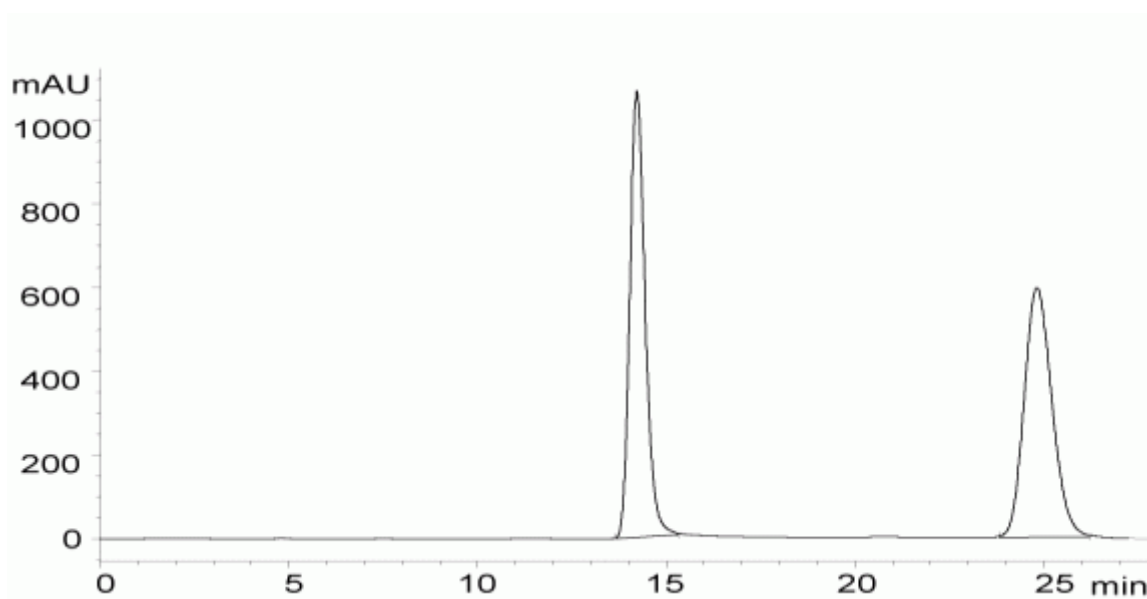

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AS-RH 150 x 4.6mm / 5µm
Mobile Phase Composition	water / acetonitrile = 60 / 40
Chromatographic Mode	Reverse Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	0.01
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 250 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	14.21	2.55	1 - 2	2.04
2	24.8	5.2		