



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Troger's Base
---------------	---------------

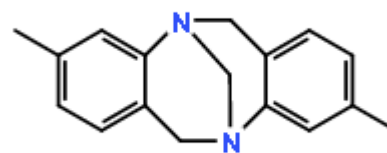

CHROMATOGRAM

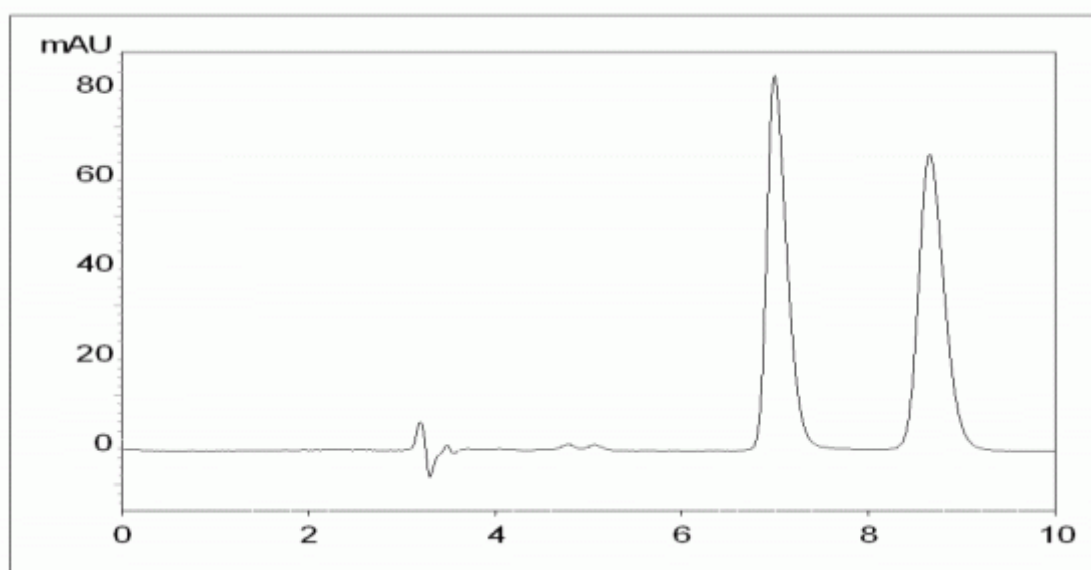

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AS-RH 150 x 4.6mm / 5µm
Mobile Phase Composition	H3PO4 aq. pH2.0 / acetonitrile = 65 / 35
Chromatographic Mode	Reverse Phase

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 254 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	6.99	2.88	1 - 2	3.71
2	8.64	3.8		