


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

DL-Tyrosine

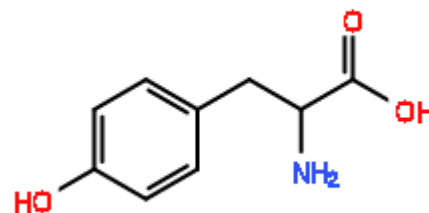

CHROMATOGRAM

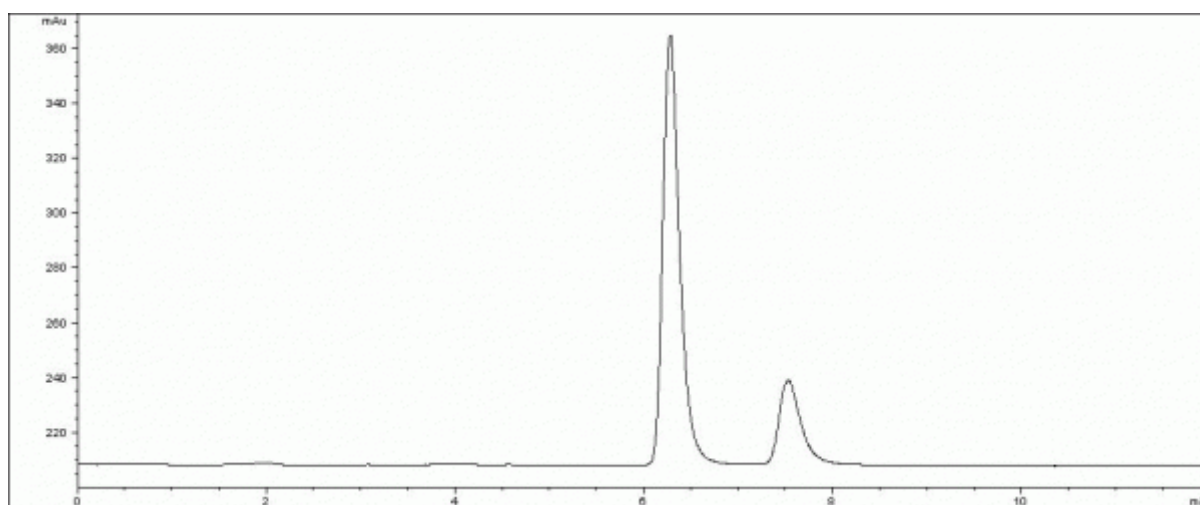

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® ZWIX(+) 150 x 3mm / 3µm
Mobile Phase Composition	25mM formic acid + 25mM diethylamine in methanol / acetonitrile / water = 49 / 49 / 2
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	6.6
Detection	ELSD	Elution order	L/D

No.	Rt	K'	Alpha Value	Resolution
1	6.28	2.6	1 - 2	3.51
2	7.53	3.31		