


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Metanephrine

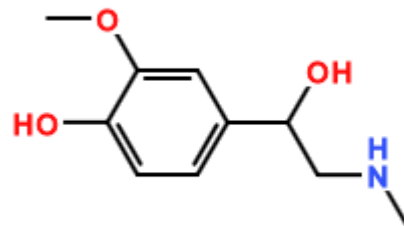

CHROMATOGRAM

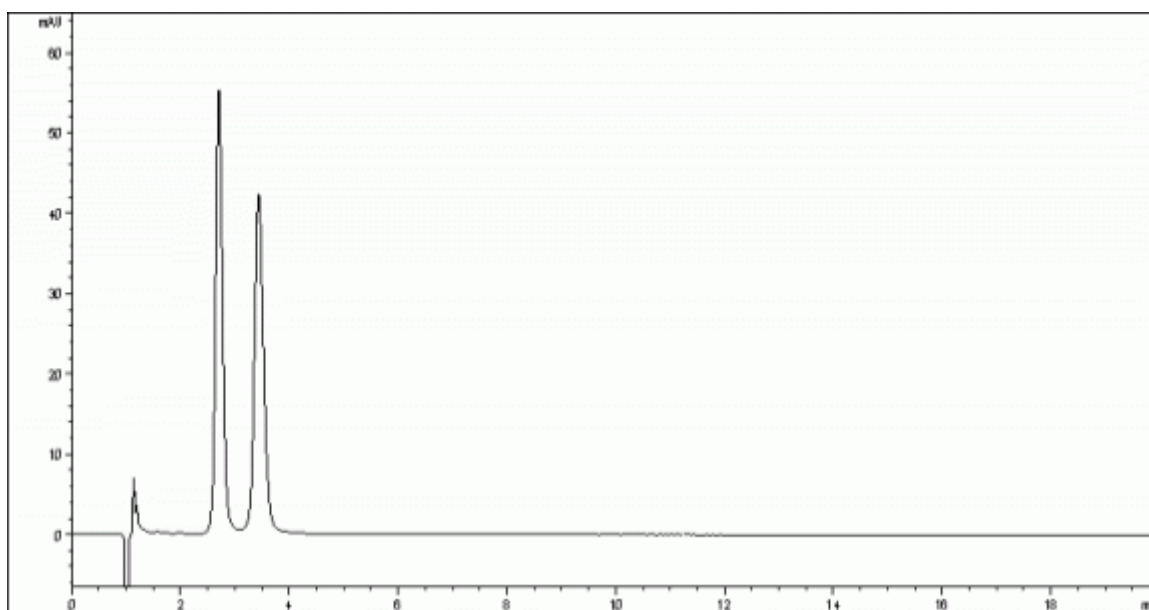

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® CBH 100 x 4mm / 5µm
Mobile Phase Composition	10mM NH4OA? aq. pH5.8 / methanol = 95 / 5
Chromatographic Mode	Reverse Phase

Flow rate (ml/min)	0.9	Injection amount (mg)	0.0005
Temperature (°C)	20	Pressure (MPa)	
Detection	UV 225 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	2.71	1.54	1 - 2	2.81
2	3.44	2.23		