



## APPLICATION RESULTS

## COMPOUND IDENTIFICATION

Compound Name	Naproxen
---------------	----------

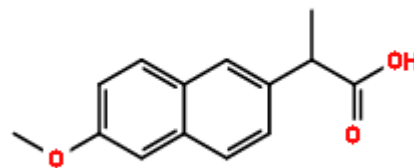

## CHROMATOGRAM

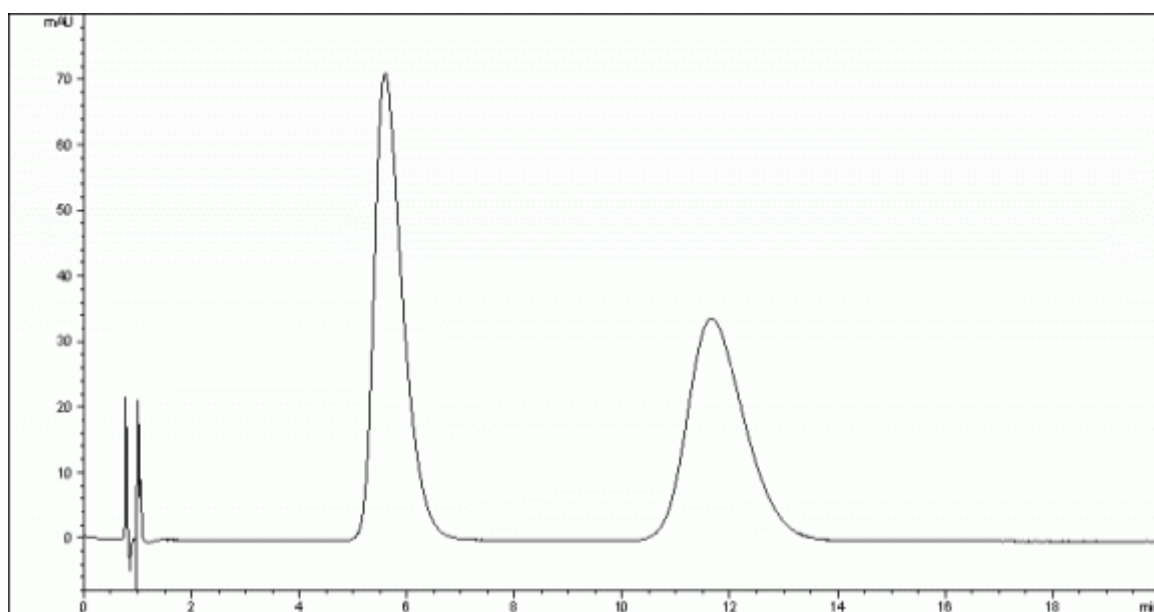

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AGP 100 x 4mm / 5µm
Mobile Phase Composition	10mM NH <sub>4</sub> OA aq. pH7.0 = 100
Chromatographic Mode	Reverse Phase

Flow rate (ml/min)	0.9	Injection amount (mg)	0.0005
Temperature (°C)	20	Pressure (MPa)	
Detection	UV 225 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	5.6	4.52	1 - 2 2.32	4.19
2	11.66	10.49		