


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Verapamil

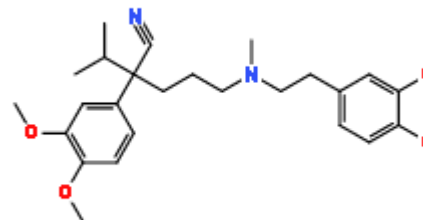

CHROMATOGRAM

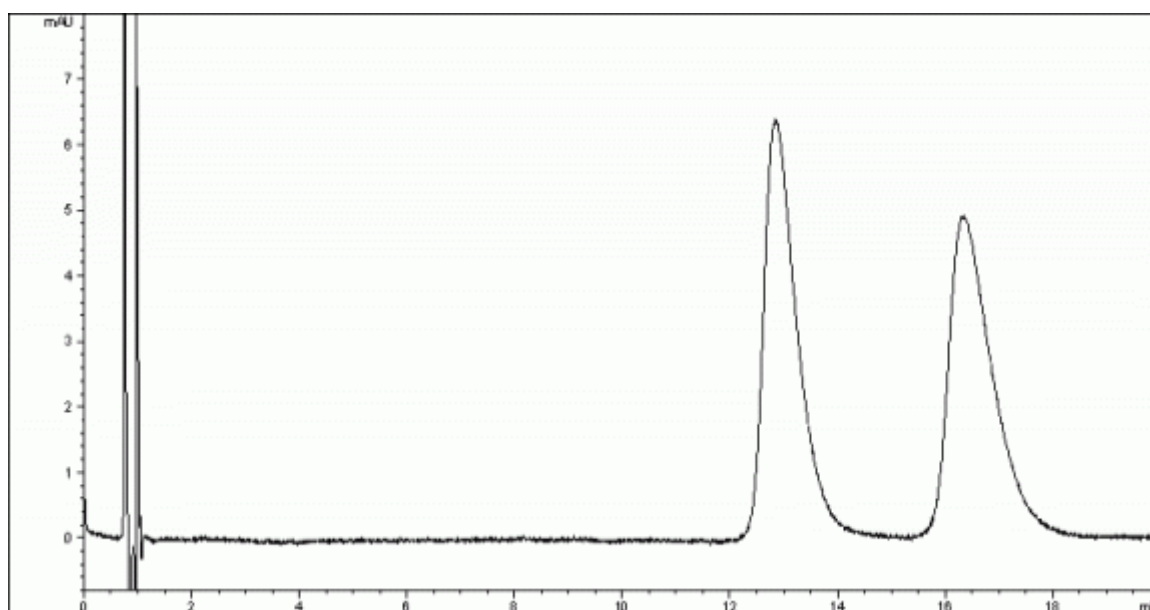

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AGP 100 x 4mm / 5µm
Mobile Phase Composition	10mM NH ₄ OA aq. pH7.0 / acetonitrile = 90 / 10
Chromatographic Mode	Reverse Phase

Flow rate (ml/min)	0.9	Injection amount (mg)	0.0005
Temperature (°C)	20	Pressure (MPa)	
Detection	UV 225 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	12.85	11.66	1 - 2	2.77
2	16.35	15.11		