



## APPLICATION RESULTS

## COMPOUND IDENTIFICATION

Compound Name	(±)-2,3-O-Benzylidene threitol
---------------	--------------------------------

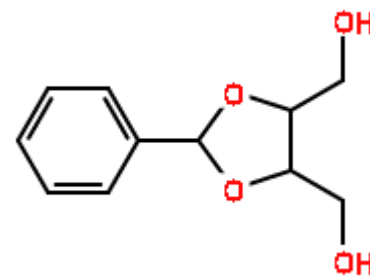

## CHROMATOGRAM

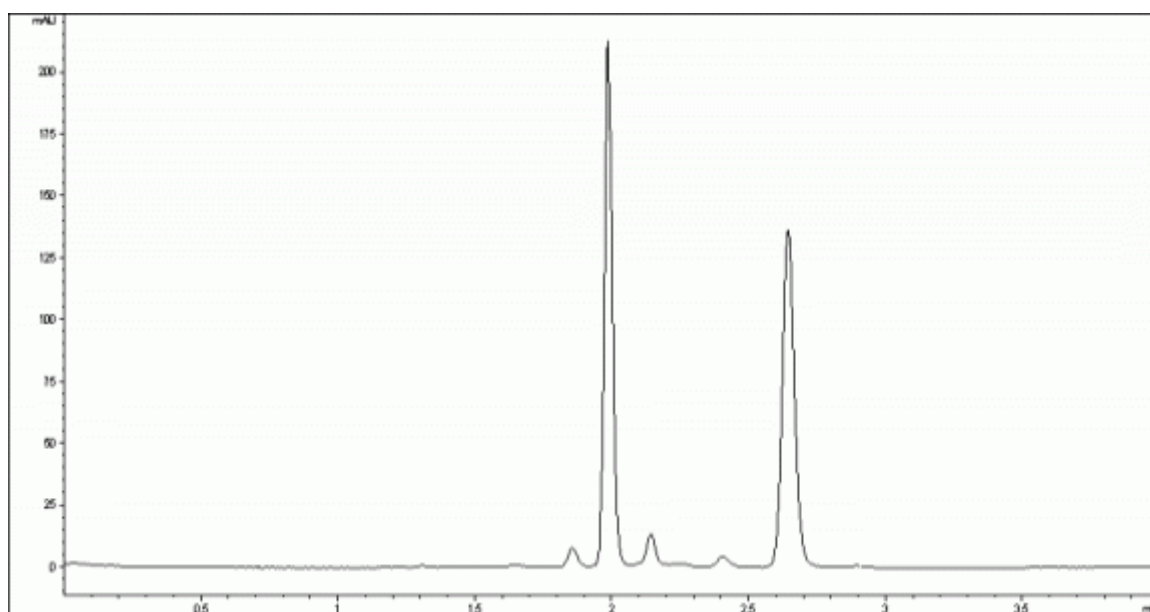

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IA-3 150 x 4.6mm / 3µm
Mobile Phase Composition	methanol = 100
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	1	Injection amount (mg)	0.001
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 210 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	1.99	0.17	1 - 2 3.33	10
2	2.65	0.55		