

APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Tetramisole



CHROMATOGRAM

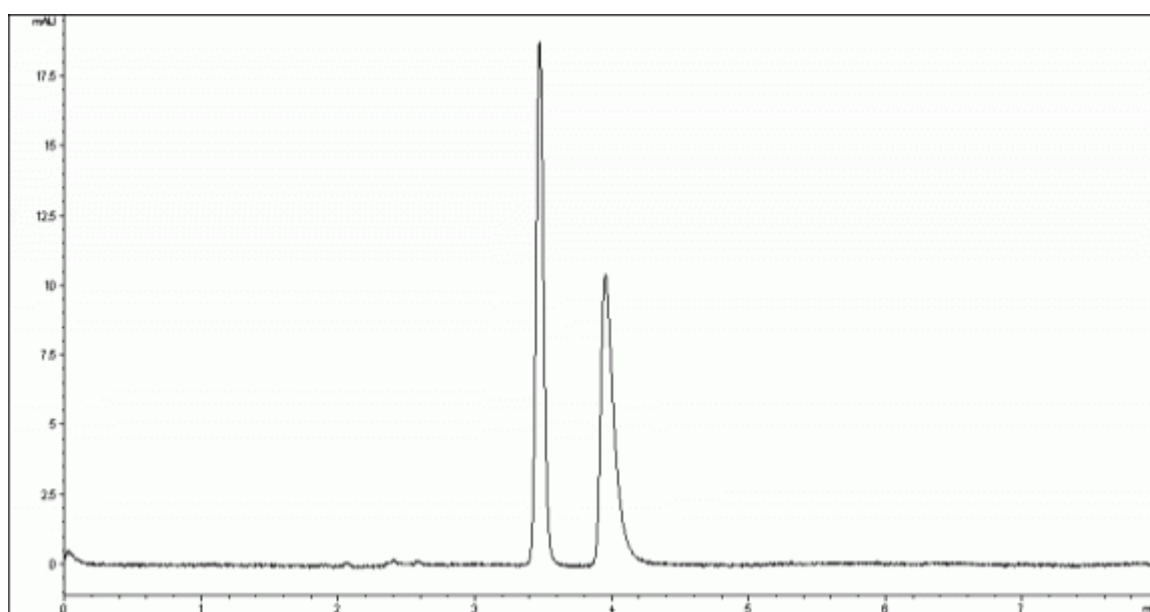

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IA-3 150 x 4.6mm / 3µm
Mobile Phase Composition	acetonitrile / ethanolamine = 100 / 0.1
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	1	Injection amount (mg)	0.001
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 254 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	3.47	1.04	1 - 2	3.61
2	3.95	1.32		