



## APPLICATION RESULTS

## COMPOUND IDENTIFICATION

Compound Name	Wieland Miescher Ketone
---------------	-------------------------

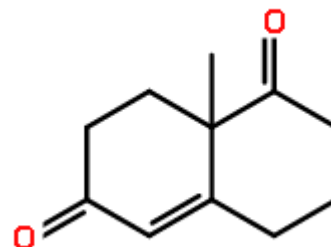

## CHROMATOGRAM

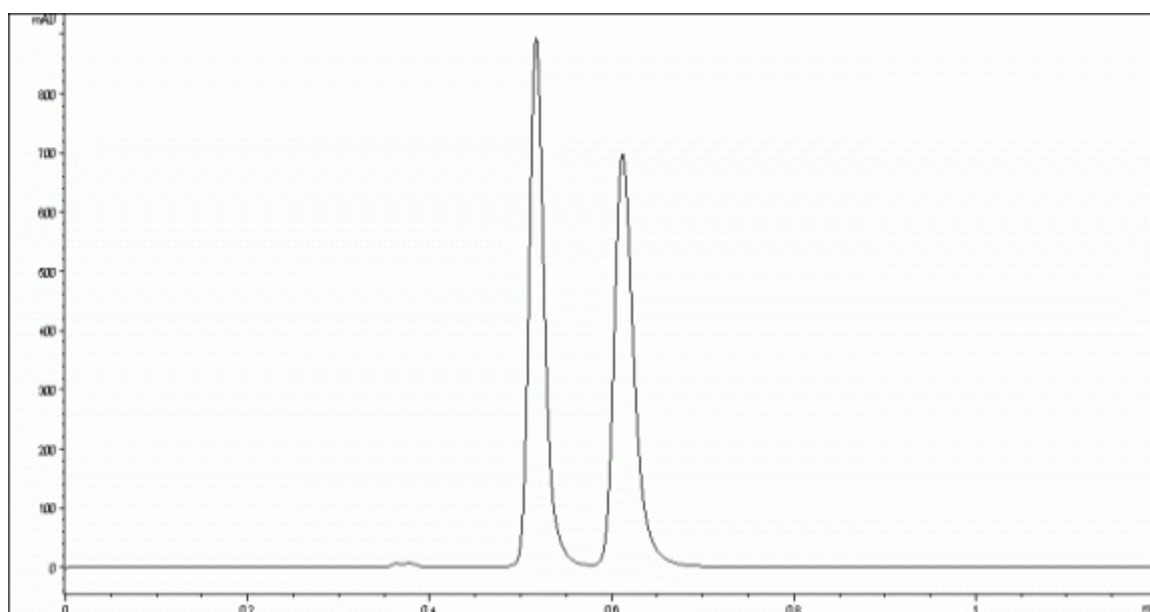

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IA-3 50 x 4.6mm / 3µm
Mobile Phase Composition	methanol = 100
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	2	Injection amount (mg)	0.002
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 254 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	0.52	0.8	1 - 2	2.84
2	0.61	1.13		