



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	4-Benzoyloxy-2-azetidinone
---------------	----------------------------

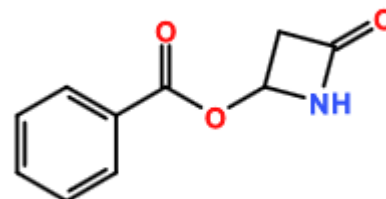

CHROMATOGRAM

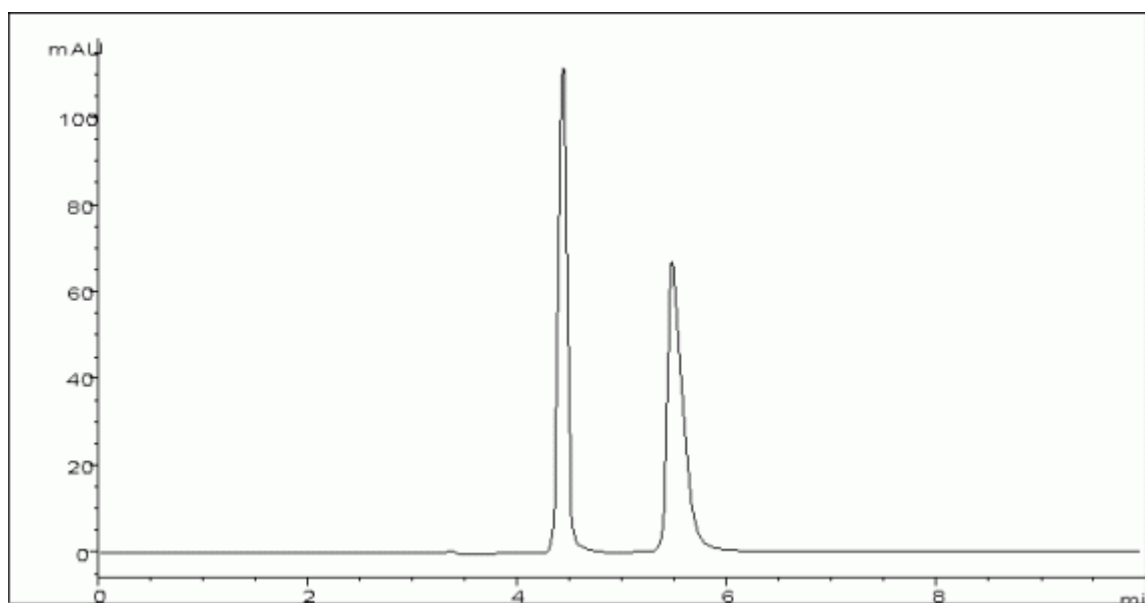

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IA 250 x 4.6mm / 5µm		
Mobile Phase Composition	acetonitrile = 100		
Chromatographic Mode	Polar Phase		

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 260 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	4.43	0.51	1 - 2	5.39
2	5.48	0.86		