


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Mirtazapine

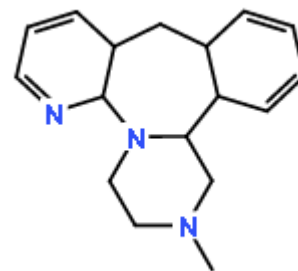

CHROMATOGRAM

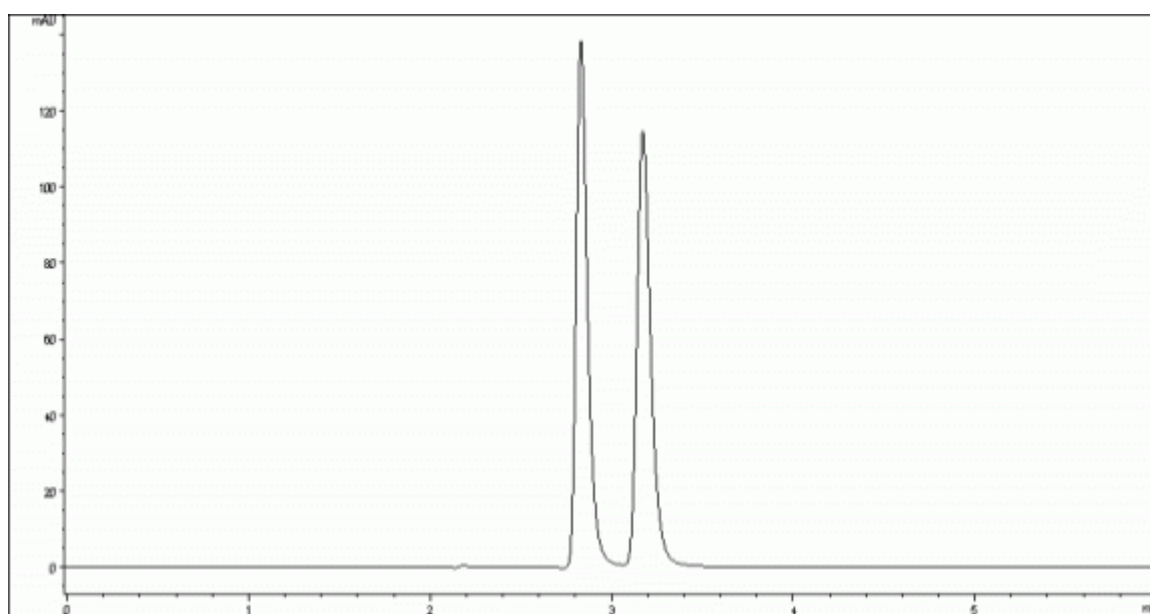

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IC 150 x 4.6mm / 5µm
Mobile Phase Composition	acetonitrile / 2-propanol / diethylamine = 95 / 5 / 0.1
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 270 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	2.83	0.48	1 - 2	2.81
2	3.17	0.66		