


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Orphenadrine

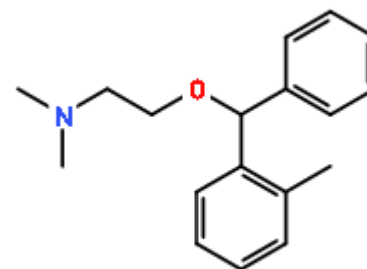

CHROMATOGRAM

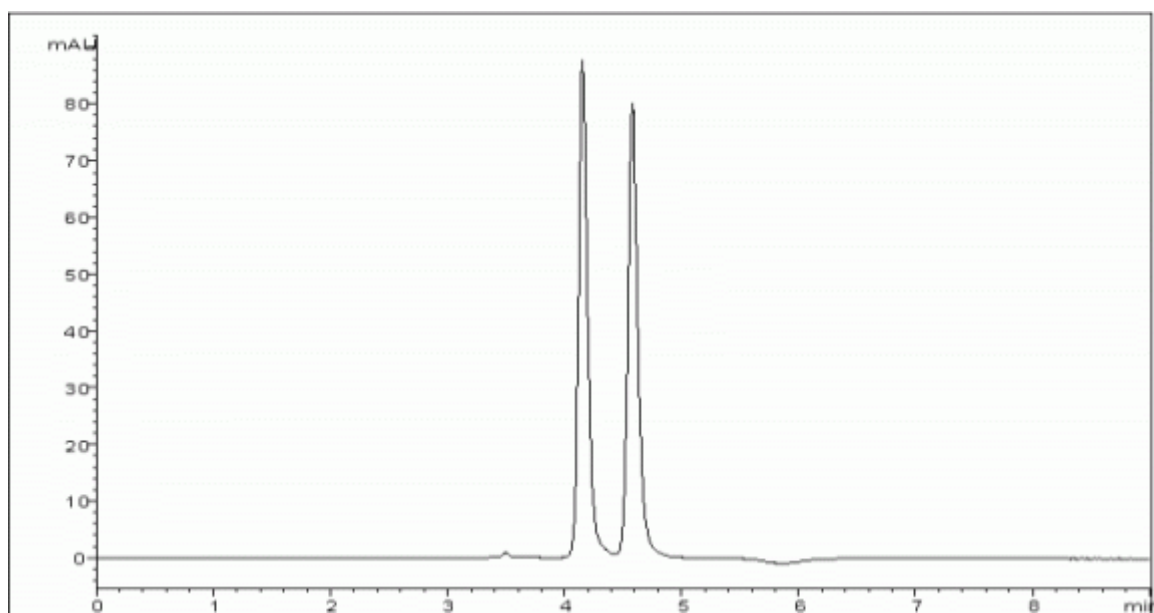

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IB 250 x 4.6mm / 5µm
Mobile Phase Composition	acetonitrile / diethylamine = 100 / 0.1
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	1	Injection amount (mg)	0.01
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 254 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	4.16	0.41	1 - 2	3.06
2	4.58	0.55		