



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Vildagliptin (racemic)
---------------	------------------------

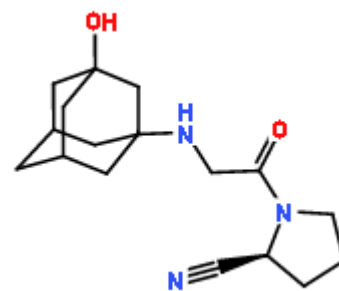

CHROMATOGRAM

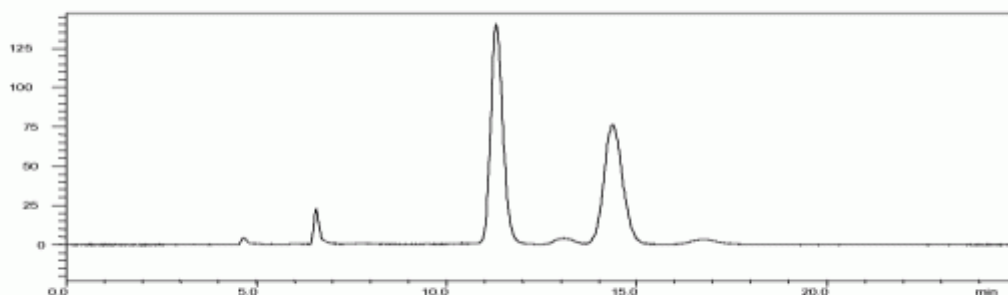

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IC 250 x 4.6mm / 5µm
Mobile Phase Composition	ethanol / diethylamine = 100 / 0.1
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 210 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	11.32	0.89	1 - 2	4.02
2	14.34	1.39		