



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Kynurenine
---------------	------------

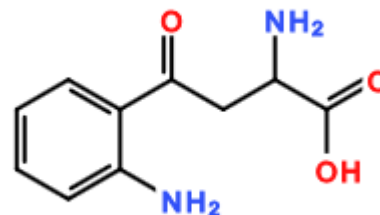

CHROMATOGRAM

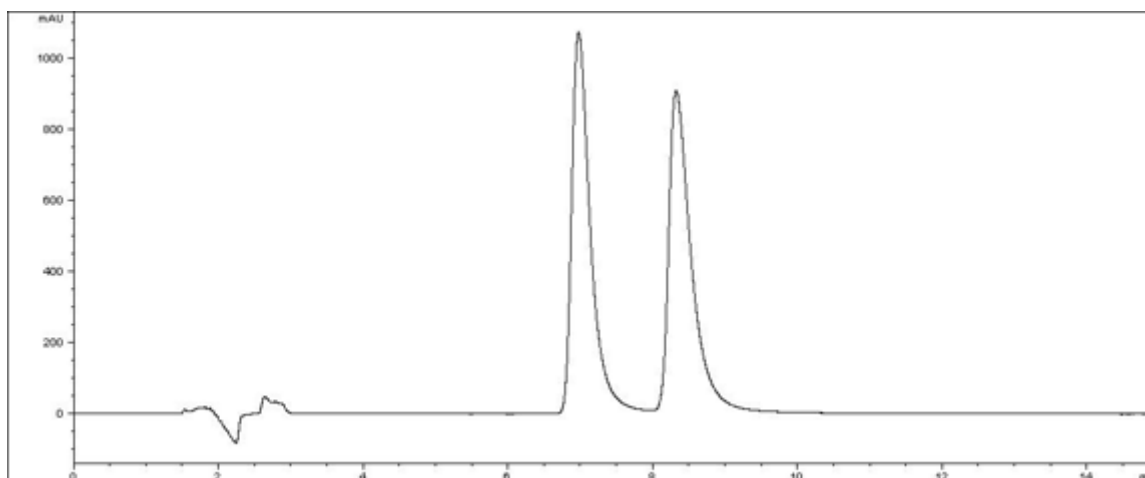

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® ZWIX(+) 150 x 3mm / 3µm
Mobile Phase Composition	25mM formic acid + 25mM diethylamine in methanol / acetonitrile / water = 49 / 49 / 2
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 230 nm	Elution order	L/D

No.	Rt	K'	Alpha Value	Resolution
1	6.97	3.38	1 - 2	2.79
2	8.32	4.23		